



STREAMLITE **LEAD**

Protection You Deserve

MADE IN THE USA



800-674-7831

READY TO PROTECT

RTP SERIES

Available in Frontal Apron and Two Piece



ANTIMICROBIAL FABRICS

R-Series aprons are made from our premium, easy-to-clean, antimicrobial Satinguard performance fabric.



READY TO SHIP IN 24 HOURS

Skip the lead times! Your R-Series apron will be ready to ship in as little as 24 hours.



NO ORDER FORMS REQUIRED

No complicated paperwork here. Just more free time.



CUSTOM NAME TAG

Want to personalize your apron? A clear, vinyl sleeve on the apron shoulder holds an embroidered name tag.



HEIGHT	5'4	5'8	5'11	6'2
CHEST	42"	45"	48"	53"
WAIST	34"	38"	42"	53"
HIPS	36"	40"	44"	48"
FRONT LENGTH	34"	37"	40.5"	42.5"
FRONT WIDTH	21"	22"	24"	26"
VEST LENGTH	22"	24"	25"	25"
SKIRT LENGTH	20"	20"	21"	21"
SIZE	SM	MED	LG	XL



FLEXBACK APRON

0.50mm LE Protection in the Front

Light Lead, Velcro Thyroid Collar at Additional Cost

Available in STREAM-LITE

Features: Elastic Back, Extra Shoulder Padding, Pocket

Sizes: S, M, L, XL

VEST SKIRT APRON

0.50mm LE Protection in the Front

Light Lead, Velcro Thyroid Collar at Additional Cost

Available in STREAM-LITE

Features: Buckle Sets, Extra Shoulder Padding, Pocket

Sizes: S, M, L, XL

QUICK SHIP


T-SERIES SEMI CUSTOM SHIPS IN 3-4 WEEKS

Our Quick Fit measuring skin provides easy sizing for your convenience. **NO TAPE MEASURES!**

Available in sizes Petite to 2XL and .35MM or .5MM protection.

Available in all two weight options for your convenience XSTREAMLITE or STREAMLITE coverage.

***CHOOSE FROM SEMI CUSTOM COLOR CHART**

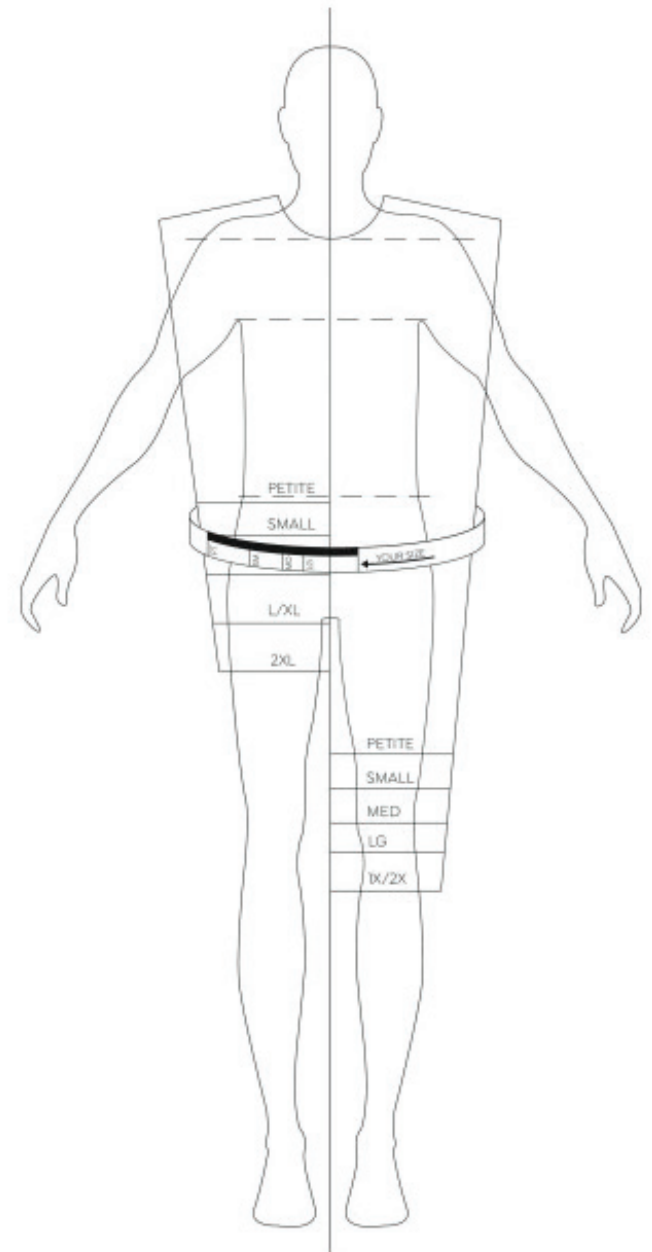

Streamlite's EZ Comfort-Flex Front Aprons offer either a buckle or Velcro®-style hook & loop closures, a comfortable unisex fit, dozens of fabric options, and free embroidery!

*CHOOSE FROM THE T-SERIES COLOR CHART



Optional Add-On: Thyroid Shield

*with binding or without



Streamlite's EZ Full Wrap Vests feature a secure Velcro®-style overlapping hook & loop closure, a comfortable unisex fit.

The EZ Full Wrap Skirts feature a combination buckle and Velcro®-style hook & loop secure closure, a comfortable unisex fit and free embroidery!

CUSTOM

FRONT PROTECTION



*CHOOSE SIZE + COLORS FROM CUSTOM CHART

P-FLEX BACK APRON

Details

Representing the gold standard of quality, the Flexback apron design provides enhanced back relief, comfort, and excellent weight distribution.

Criss-cross back straps are specially designed with elastic to “flex” the apron’s weight burden off your shoulders to your hips, thereby providing support to your lower back.

A velcro backstrap (shown) is included and secures the apron around your shoulder area to prevent sliding.

XSTREAMLITE	LEAD FREE / LIGHTWEIGHT
STREAMLITE	STANDARD LEAD FREE

FRONTAL APRON FEATURES

- **Elastic Back**
Provides enhanced adjustability and lumbar support. Shifts weight from shoulders & hips.
 - ▶ **Shoulder Padding**
Premium foam padding in the shoulders provides extra comfort.
 - ▶ **Apron Pocket**
 - ▶ **Velcro Closure**
- **Backstrap**
Velcro backstrap securely fastens shoulder area.

ADDITIONAL CUSTOMIZATION OPTIONS

- ▶ THYROID COLLAR
- ▶ EXTRA POCKET
- ▶ LOGO
- ▶ DROP SHOULDER
- ▶ EMBROIDERY
- ▶ SLEEVE
- ▶ WINGS
- ▶ BACK MESH

P-SURGIDROP

Maintain your sterile field with the quick release Surgical Drop Apron. Easy-release velcro closures allow the surgical drop to fall even when underneath your scrubs. With padded shoulders and a criss-cross back strap design the surgical drop apron provides maximum comfort for all users. Hanging straps are included to allow for easy storage on any apron rack.

ADDITIONAL CUSTOMIZATION OPTIONS

- ▶ THYROID COLLAR
- ▶ EXTRA POCKET
- ▶ EMBROIDERY
- ▶ LOGO

*** CHOOSE COLORS FROM CUSTOM CHART**





BACK RELIEF APRON

Details

Save yourself from back pain with Streamlite's Back Relief apron. This apron features a unique padded belt system constructed with a bonded elasticized panel that will not lose its shape. When engaged, the special back relief belt provides lumbar support while transferring the weight of the apron from your shoulders to your hips for increased comfort and back support throughout the life of the apron. A velcro backstrap is also included to secure the apron around your shoulder area and prevent sliding.

ADDITIONAL CUSTOMIZATION OPTIONS

- ▶ THYROID COLLAR
- ▶ EXTRA POCKET
- ▶ EMBROIDERY
- ▶ LOGO

***CHOOSE COLORS FROM P CHART**

FULL PROTECTION P-WRAP



***CHOOSE COLORS FROM CUSTOM CHART**



FEATURES

Offering full protection, this back relieving design transfers the weight from your shoulders to your hips with the aid of a 6" lumbar support belt. The Wrap Around Back Relief apron features velcro and buckle closures to ensure a proper fit. The Balance configuration features seam-to-seam, single LE protection with standard overlapping panels (see below for lead configuration).

Chest pocket included and Shoulder pads provide additional comfort. Belt is also matched to the color or pattern of your apron.

***CHOOSE COLORS FROM P-CHART**

CUSTOM

TWO PIECE PROTECTION



***CHOOSE COLORS FROM CUSTOM CHART**

TWO PIECE VEST+SKIRT

Details

MULTIPLE CHOICES IN VEST + SKIRT OPTIONS

- ▶ STANDARD VEST BUCKLE CLOSURE
- ▶ REVERSE VEST
- ▶ FLEX BACK VEST

- ▶ STANDARD SKIRT WITH POCKET
- ▶ STANDARD SKIRT WITH POCKET + SUSPENDERS

XSTREAMLITE	LEAD FREE / LIGHTWEIGHT
STREAMLITE	STANDARD LEAD FREE

VEST + SKIRT FEATURES

STANDARD VEST:

Double buckle with hook and loop closer to stay secure

- ▶ Shoulder Padding
Premium foam padding in the shoulders provides extra comfort.

REVERSE VEST:

Elastic Back

Provides enhanced adjustability and lumbar support. Shifts weight from shoulders & hips.

- ▶ Apron Pocket
- ▶ Velcro Closure

FLEX BACK VEST:

Flexback shifts 70% of the weight directly to the hips with added support of an adjustable, elastic back relief belt.

ADDITIONAL CUSTOMIZATION OPTIONS

- ▶ THYROID COLLAR
- ▶ EXTRA POCKET
- ▶ LOGO
- ▶ DROP SHOULDER
- ▶ EMBROIDERY
- ▶ SLEEVE
- ▶ WINGS
- ▶ BACK MESH

APRON STRAPS FOR STORAGE



NON-SLIP FABRIC FOR COMFORT





SEMI CUSTOM COLORS

STANDARD NYLON



11 - Royal Blue




12 - Navy Blue



20 - Green



21 - Teal

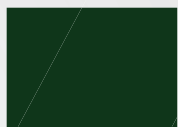

22 - Forrest



23 - Grey



33 - Black



41 - Sunshine



43 - Orange



51 - Red



52 - Hot Pink



53 - Purple




54 - Berry

COMFORT NYLON



91 - Stratos

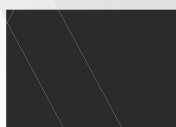

93 - Midnight

REINFORCED NYLON



61 - Sapphire



62 - Ruby



63 - Emerald




64 - Onyx

SPECIALTY NYLON




69 - Skulls

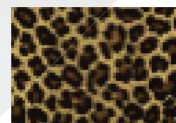


71 - Leopard




72 - Lots O Dots

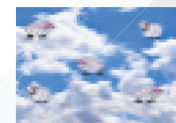

73 - Sleepy Sheep

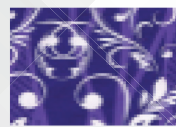

77 - Life Vines

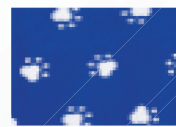


76 - Blue Paws




78 - Waterfall




80 - Red Paws

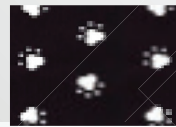


81 - Black Paws

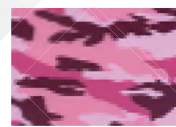

82 - Pink Camo

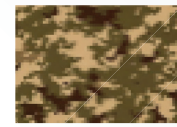

85 - Digital Camo

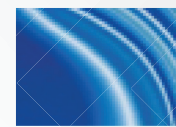

92 - Ocean

BINDING COLORS



K - Black



G - Gray


ALL FABRICS ARE LIGHTWEIGHT, LIQUID REPELLENT AND FILAMENT REINFORCED
***CHOOSE TWO OPTIONS FOR INSIDE AND OUTSIDE / FIRST CHOICE & SECOND CHOICE**

CUSTOM COLOR CHART


STANDARD NYLON + BINDING




ANTIMICROBIAL + WIPE AWAY

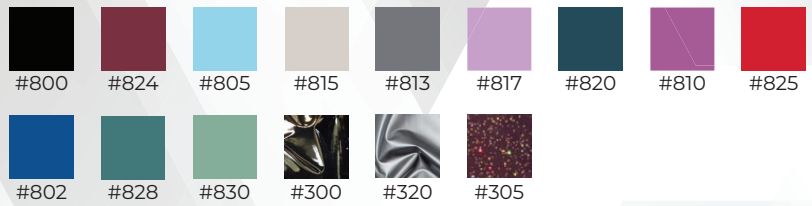


DIAMOND-REINFORCED

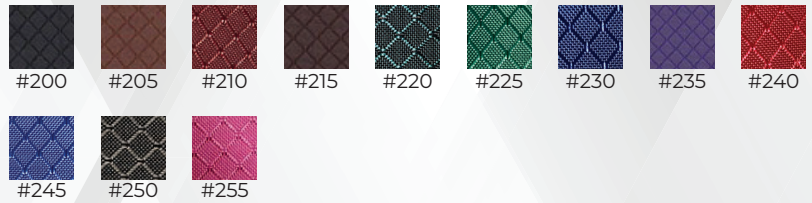

FUN FASHION NYLON

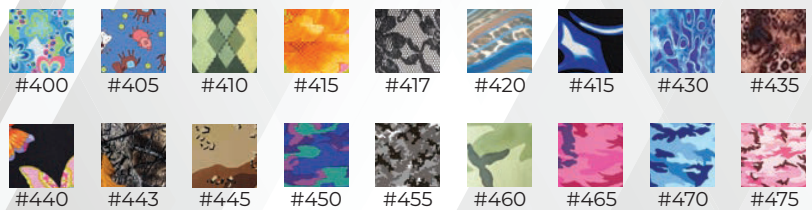

FUN FASHION NYLON

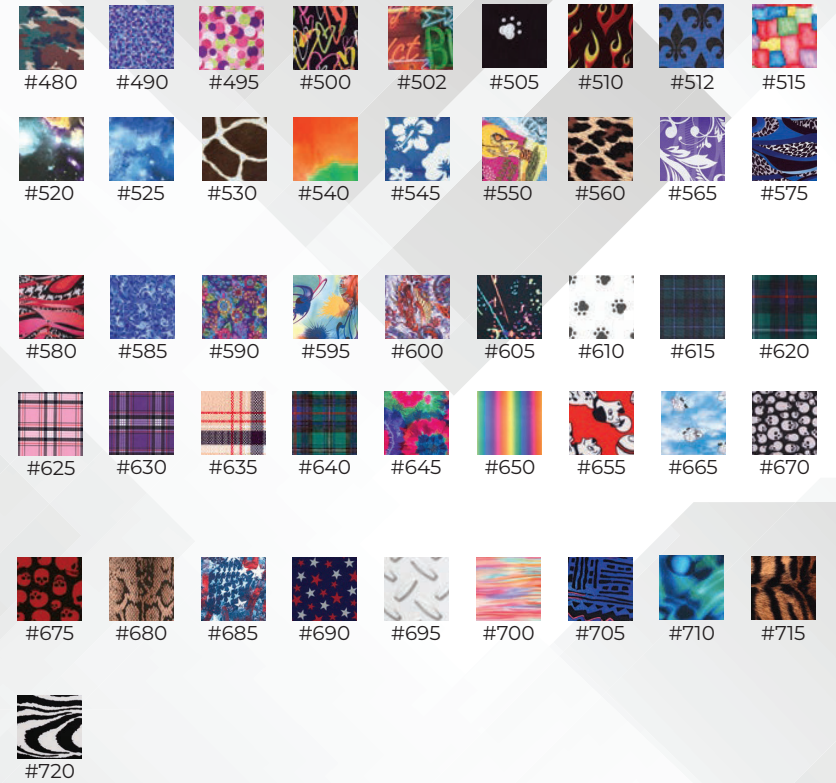


CUSTOM MEASURE GUIDE

APRON SIZE CHART

UNISEX

Small	Med	Large	XL	
5'4"	5'8"	5'11"	6'2"+	Height
42"	45"	48"	53"	Chest
34"	38"	42"	46"	Waist
36"	40"	44"	48"	Hips
34"	37"	40.5"	42.5"	Front Length
21"	22"	24"	26"	Front Width
22"	24"	25"	25"	Vest Length
20"	20"	21"	21"	Skirt Length
36"	38"	40.5"	42.5"	WABR Length

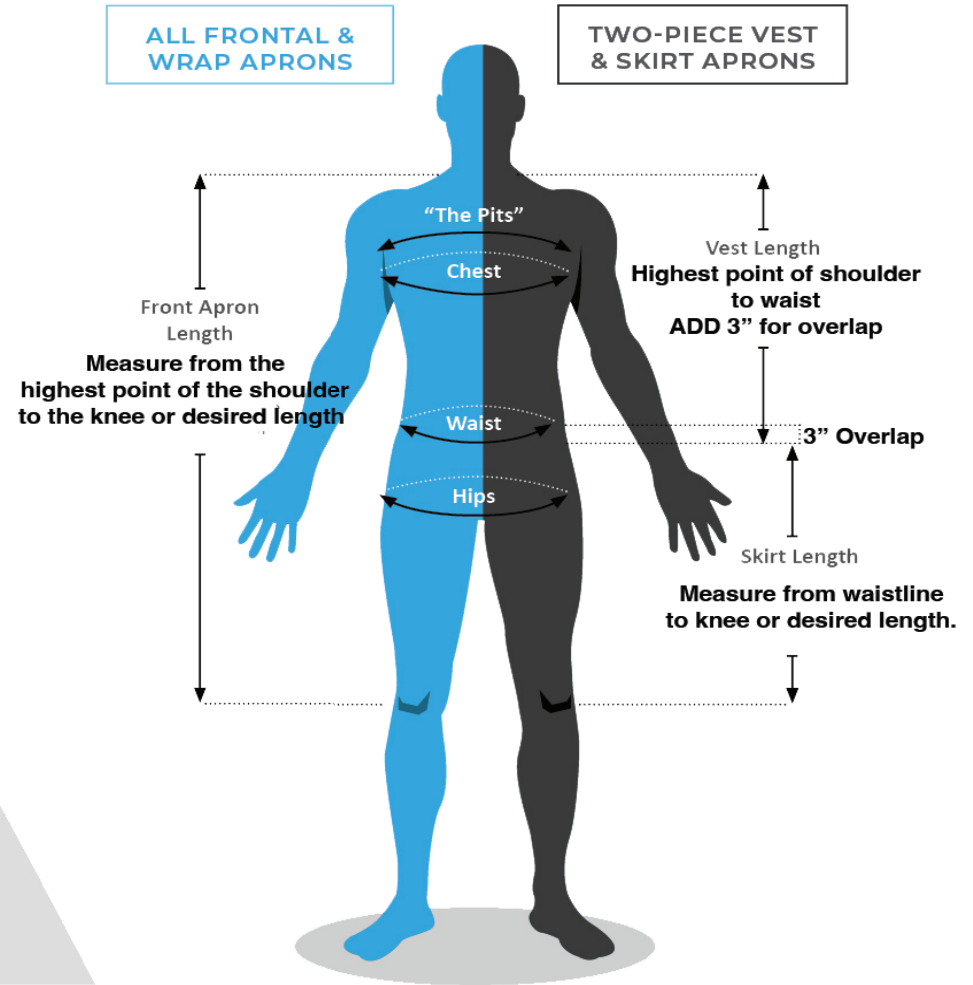

CUSTOMIZATION OPTIONS & INSPIRATION

When it comes to designing your perfect custom apron, Streamlite is #1. Our great selection of add-ons for your apron can improve your protection and comfort. Each customization is clearly defined to ensure you make the right choices for you.



APRON SUSPENDERS



SURGICAL DROP SHOULDERS

Specialty velcro tabs can be added to different apron styles allowing for easy removal of your apron to maintain sterile field.





AXILLARY WINGS

Axillary wings are a protective flap sewn-in or “snapped-on” to your apron and offer vital protection to vulnerable breast tissue and red bone marrow areas.

SIDE EXTENSION

Side Extensions can provide an additional +2 inches of protection on both sides and provide important lateral protection.



Protection levels will match your apron specifications.



APRON SLEEVE

Apron sleeves wrap around your entire arm and offer comprehensive protection of vulnerable breast tissue and red bone marrow areas. Sleeves come standard in 0.25mm LE and can be sewn-in permanently or offered in our removable “Snap-on” option.

EXTRA POCKETS

Need an extra pocket or two? Whether you need one on your skirt, inside your vest or even on your knees (yes, we did it), we can accommodate most requests. Shown to the right is a spacious L-shaped pocket which can be added to your skirt.



CUSTOMIZATION OPTIONS CONTINUED

LOGO + EMBROIDERY

Personalize your apron with your name or logo anywhere on your apron. Common embroidery locations are on the apron pocket and on a matching thyroid.



VELCRO PATCH

A customized, embroidered velcro patch can be added to a specific area on your apron. Typical locations are on the pocket, above pocket, or to the left of your pocket.



LUXTREND FABRIC

Have your apron tailored to perfection with LuxTrend's patterned accents on the apron sides, pocket and thyroid collar.

If you have a special fabric you'd like to use, we can use that too!



NON-SLIP FABRIC

A high friction, anti slip fabric that is added to the back interior of your skirt or other key areas in different apron styles which reduces slippage and improves fit.



APRON STORAGE

EXTEND THE LIFE OF YOUR APRON



WALL RACK

**CAPACITY: 5 APRONS
PER CONSOLE**

Prolong the life of your radiation protective garment with our apron wall rack console. Each console holds up to 5 protective aprons with the combined left and right consoles storing up to 10 aprons. The wall console features hooks on the top of the hangers for lead skirts. *Does not include mounting hardware.



PEG RACK

CAPACITY: 1-9 APRONS

The apron Peg rack is made of steel construction and easily screws into wall studs. It can store up to 9 aprons but comes in a variety of sizes that accommodate the storage of 1, 2, 4, 7, or 9 lead aprons. This affordable wall storage for aprons is essential to prolong the life of your protective aprons. *Mounting hardware is not included.

When not in use, garments should be placed on a hanger or stored using Streamlite's apron storage equipment. Avoid contact with hot surfaces or exposure to direct sunlight. Also avoid contact with sharp objects, as these may damage the garment's protective core material.

For transportation, garments should be carefully rolled up, not folded, and placed into a box or Streamlite duffel bag.

PORTA RACK



CAPACITY: 8 APRONS



Ideal for use in the O.R., Specials, and Cath Labs, this portable apron storage unit is great quality and affordable. Our mobile lead apron rack has a 5-leg, caster-equipped base, making it the perfect portable apron rack to store and transport up to 8 aprons. This solid aluminum rack is sturdy and stable, and can be quickly moved out of the way.

APRON STORAGE

SMALL & COMPACT APRON STORAGE

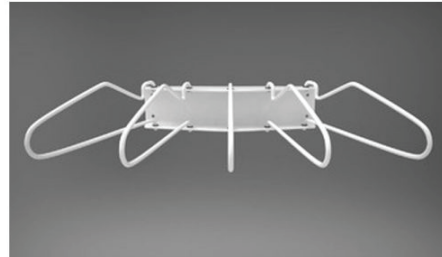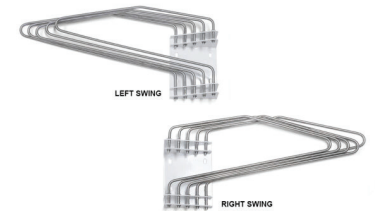


APRON STORAGE

MOBILE RACK OPTIONS

MIGHTY RACK

CAPACITY: 10 APRONS

The Mighty-Rack mobile apron rack's solid steel design holds up to 10 aprons. The pedestal base has rolling wheels for easy portability. Convenience hooks on each hanger can hold lead skirts.

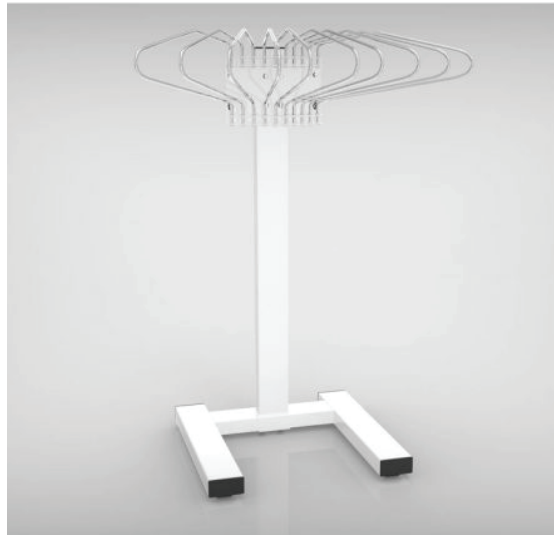

GARMENT RACK

CAPACITY: 10 ARM RACK

The Deca-Rak has 10 pivoted arms that fan out to 180 degrees. It is mounted on locking casters and has a high-tech finish of chromium and white that enhance any room decor. Specifications: 54" H x 52" W x 24" D With center arm at 90°, angle from racks, depth is 26" Weight Capacity: approximately 600 lbs.



DUFFLE BAG



APRON STORAGE

MOBILE RACK OPTIONS CONT.

VALET RACK

CAPACITY: 10 APRONS

The Valet Rack mobile apron rack's comes equipped with five closed hangers and holds up to 300 lbs. It is sturdy yet easily transported between work spaces. This apron rack is equipped with heavy-duty locking casters and is extremely mobile allowing aprons, thyroid shields, etc. to be transported easily to workstations.

Constructed of steel and finished with a white powder coating. Five end-flared closed hangers are included with each unit with the option to add up to five more.


1 year warranty against manufacturing defects.

Available in 36" and 48" wide.




INSET SHOWN:

Unit overall size: 58" High x 24" Deep x 48" Wide

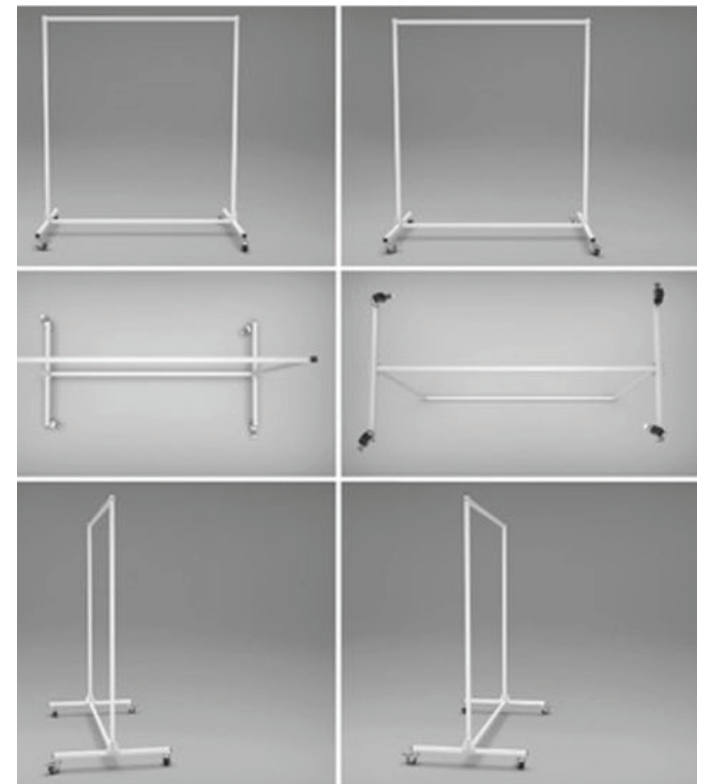

GLOVES

PERSONAL RADIATION PROTECTION

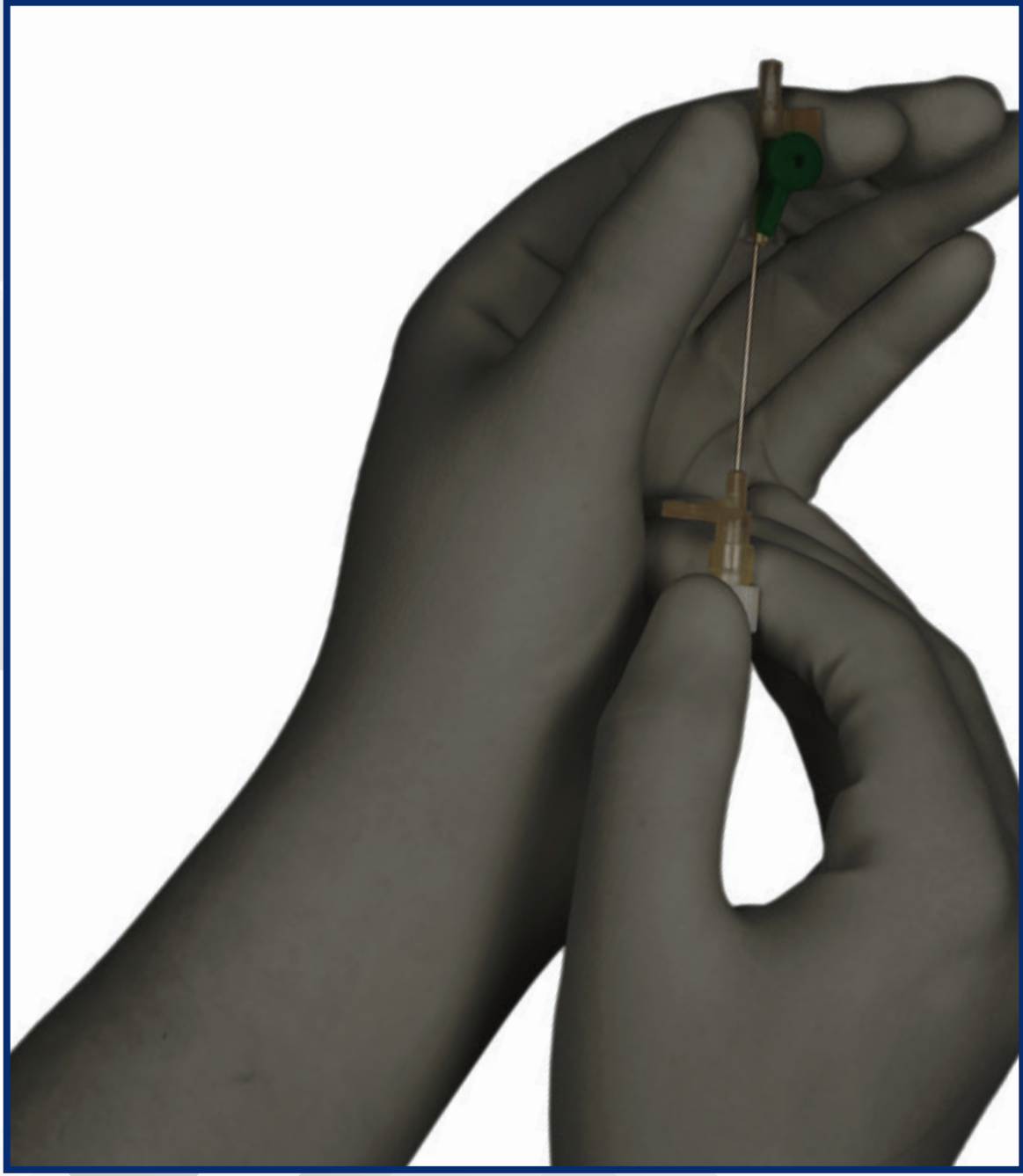

"By simply holding the head of a trauma patient four times per week with unshielded hands during cervical spine radiography, a healthcare worker could receive more than 1000 mSv of radiation (that's twice the annual limit of radiation exposure). Over-exposure to radiation can have debilitating repercussions such as cancer, and, let's be honest, nobody wants that. Say no to cancer and protect yourself with radiation attenuating gloves.

ATTENUATION GLOVES


OPTIONS



Latex Free + Lead



Latex Free + Lead Free

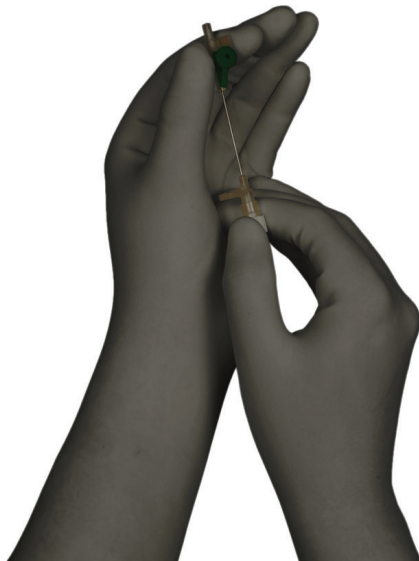

AVAILABLE IN LEADED AND NON-LEADED STYLES
SIZES 6-9 WITH HALF SIZES AVAILABLE



ANGIOGRAPHIC LEAD GLOVE



ANGIO MITT



RR SLEEVE



SLIT MITT

EYEWEAR

QUALITY. SERVICE. INTEGRITY.



*"The typical interventional cardiologist gets exposed to 20 to 30 mSv each year. This level of exposure is nearly **10 times greater** than the level of exposure received by the entire body."
-- US National Library of Medicine*

WHO USES PROTECTIVE GLASSES?


Your eyes are 4 times more sensitive to radiation than other parts of your body. Not too long ago, few doctors used leaded eyewear. With the increase in occurrences of early onset cataracts, lead glasses are now a staple in fluoroscopic procedures. Protective leaded eyewear is used by roughly 95% of physicians and staff in the following vocations (but not limited to):



CATH LABS



EP LABS

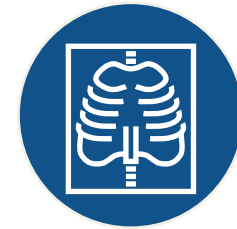

INTERVENTIONAL
RADIOLOGY



IMAGING CENTERS



ORTHOPEDICS



UROLOGY

PROTECTION SUMMARY

**Lead Equivalence - 0.50 or 0.75mm Pb X-ray Tube Voltage - 60 kV – 150 kV
IEC 61331-3:2014 Regulation (EU) 2016/425 Module B EU Type Examination**

EYEWEAR BRANDS



We carry many brands and styles (including pediatric). Enclosed are a few samples of the top sellers. If you wish to see the full line we offer, please reach out to your local representative for additional options!





BARRIERS

LEAD BARRIERS AND SHIELDS



MOBILE BARRIERS

We represent a comprehensive line of x-ray protective barriers, windows, and other leaded shield products. From fully custom lead doors to unique windows and overhead barriers, we have the expertise to meet your needs. Many of the mobile barriers on the market restrict viewing due to large opaque panels. The line we represent has pioneered a new method of producing barriers. First, instead of glass, they use a light-weight lead acrylic which is substantially more durable. The acrylic makes up most of the barrier with only a small opaque panel. This new, clean looking design is lightweight and provides an excellent field of vision. In line with our commitment to quality, most barriers and components are produced in the USA.

FULL MOBILE BARRIER



HALF MOBILE BARRIER



Diagnostic x-ray procedures frequently require working personnel to remain in the examination room where they are subject to exposure from secondary radiation.

The clear leaded acrylic barrier can provide total-body x-ray shielding during extended procedures. Our clear acrylic barrier provides 0.50mm LE protection and rests on an opaque lower panel with an equivalent 0.50mm LE.

* AVAILABLE IN FULL AND HALF ACRYLIC VIEWS

Smooth-rolling casters give the barrier great mobility. The acrylic lead barrier is shipped unassembled but can be easily assembled using a screwdriver and wrench.

**The base features 0.50mm Pb lead sheeting
Acrylic and Aluminum/Steel are made in the USA**

MOBILE LEADED BARRIERS

VARIOUS STYLES & SIZES AVAILABLE



SHS-LB1224

2.0mm Pb
Lead Glass window: 12" x 24"
Overall dimensions: 74" H x 52½"
W X 25" D



SHS-LB2024-N

Overall size: 75 (h) x 27.5 in (w)
Window panel: 20 (h) x 24 in (w)
Opaque panel: 47 (h) x 27.5 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB2436

2.0mm Lead equivalent opaque
panels of solid leaded x-ray
protection with contoured sides.
24" x 36" window for easy viewing.
Overall size: 74" H x 39" W x 25" D

*CHOICE OF GLASS OR ACRYLIC



SHS-LB2430-MRI

MRI-Safe Mobile Leaded
Aluminum Barrier with Window
Overall size: 75 (h) x 33.5 in (w)
Window panel: 24 (h) x 30 in (w)
Opaque panel: 48 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC

***AVAILABLE INSERT OPTIONS**
2.0mm lead equivalent high quality Schott leaded glass
0.5mm lead equivalent leaded acrylic

MOBILE LEADED BARRIERS

VARIOUS STYLES & SIZES AVAILABLE



SHS-LB2430-B

Tilted Mobile Leaded Barrier with
24" Window
Overall size: 70" H x 34" W X 25" L
Window panel : 24 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB2436-S

"SHORTY" designed for shielding
nuclear medicine personnel.
Overall size: 64"H x 39"W
Window panel: 24"H x 36"W
Opaque panel: 34"H x 36"W

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3648

Overall size: 75 (h) x 52.5 in (w)
Window panel: 30 (h) x 48 in (w)
Opaque panel: 36 (h) x 48 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3648-MRI

MRI-Safe Mobile Aluminum Barrier
with 30" x 48" Window
Overall size: 75 (h) x 52.5 in (w)
Window panel: 30 (h) x 48 in (w)
Opaque panel: 36 (h) x 48 in (w)

*CHOICE OF GLASS OR ACRYLIC

***AVAILABLE INSERT OPTIONS**
2.0mm lead equivalent high quality Schott leaded glass
0.5mm lead equivalent leaded acrylic

MOBILE LEADED BARRIERS

VARIOUS STYLES & SIZES AVAILABLE



SHS-LB3672

Overall size: 75 (h) x 75.5 in (w)
Window panel: 36 (h) x 72 in (w)
Opaque panel: 36 (h) x 72 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3036

Overall size: 75 (h) x 33.5 in (w)
Window panel: 36 (h) x 30 in (w)
Opaque panel: 22 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3048

Overall size: 75 (h) x 33.5 in (w)
Window panel: 48 (h) x 30 in (w)
Opaque panel: 22 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3048-MRI

MRI-Safe Mobile Leaded Aluminum Barrier with Window
Overall size: 75 (h) x 33.5 in (w)
Window panel: 48 (h) x 30 in (w)
Opaque panel: 22 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC

***AVAILABLE INSERT OPTIONS**
2.0mm lead equivalent high quality Schott leaded glass
0.5mm lead equivalent leaded acrylic

MOBILE LEADED BARRIERS

VARIOUS STYLES & SIZES AVAILABLE



SHS-LB3060

Overall size: 75 (h) x 33.5 in (w)
Window panel: 60 (h) x 30 in (w)
Opaque panel: 12 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB3060-MRI

MRI-Safe Mobile Leaded Aluminum Barrier with Window
Overall size: 75 (h) x 33.5 in (w)
Window panel: 60 (h) x 30 in (w)
Opaque panel: 12 (h) x 30 in (w)

*CHOICE OF GLASS OR ACRYLIC



SHS-LB7429

A screen that can be easily raised and lowered by hand
0.50mm lead equivalent light weight leaded acrylic window
2.0mm lead equivalent steel and lead frame
Overall size: 70 (h) x 31 in (w)



SHS-LB8060

The barrier may be rotated vertically 360 degrees and has a flexible Velcro attached, lead rubber shield that is 14 inches wide.
0.5mm lead equivalent
Overall size: 31.5 (h) x 23.6 in (w)

*AVAILABLE INSERT OPTIONS

2.0mm lead equivalent high quality Schott leaded glass
0.5mm lead equivalent leaded acrylic

MOBILE LEADED BARRIERS

ADJUSTABLE OVERHEAD BARRIERS



SHS-PTO-001

Square Arm Overhead Lead Acrylic
Mobile Barrier with Torso Cut-Out

Window panel: 32 (h) x 24 in (w)



SHS-PTO-003

Round Arm Overhead Lead Acrylic
Mobile Barrier with Torso Cut-Out

Window panel: 32 (h) x 24 in (w)



SHS-PTO-002

Square Arm Overhead Lead Acrylic
Mobile Barrier with
0.50mm PbEq Lead Curtain

Window panel: 24 (h) x 16 in (w)
Curtain panel: 24 (h) x 14 in (w)



SHS-PTO-004

Round Arm Overhead Lead Acrylic
Mobile Barrier with
0.50mm PbEq Lead Curtain

Window panel: 24 (h) x 16 in (w)
Curtain panel: 24 (h) x 14 in (w)

PORTA SHIELDS - ADJUSTABLE SHIELDS

ADJUSTABLE BARRIER

Our adjustable mobile radiation barrier provides maximum shielding from x-ray exposure. Ideal for patient & personnel radiation protection our lead shield is easy to maneuver, and much lighter than the full & half acrylic barriers.

This enables greater mobility while maintaining a solid construction with high quality, durable wheel casters. You can also customize your barrier with your company's logo or embroidery. Velcro-lined for improved adjustability!



PORTA SHIELD

The graduated porta shield is a mobile-adjustable lead shield that provides a variety of different adjustable heights blended with a lightweight design. The "H" shaped wheels provide stability and portability.

The shield is attached by velcro which enables easy removal for cleaning.

* CUSTOM SIZING AVAILABLE



GRADUATED PORTA SHIELD

The graduated porta shield is similar to the porta shield but differentiates itself with its individual lead flaps. The flaps provide excellent versatility and usage. They are attached by velcro which enable easy removal for cleaning.



PORTA SHIELDS - ADJUSTABLE SHIELDS

VARIOUS STYLES & SIZES AVAILABLE

.....



MAX MOBILE LEAD SHIELD - H-BASE

0.5mm Pb equivalency
Adjustable height: up to 60"
Base Overall: 24" W x 24" D
H Base Dimensions: 20" W x 16" D
All steel H-base for stability. Hardened chrome steel post.

Available in two sizes:
36x36 & 48x48



STANDING MOBILE SHIELD

Compact design and opaque steel sides, 2.0mm Pb equivalency

Overall: 48" H x 21.5" W x 16" D
Easy Grip Handles
Steel H Base



SHS-LB-MB200404-PVC

HEAVY DUTY MOBILE SHIELD

Height adjustable with tightening knobs. Top of Lead Shield is 36" H when lowered and 58" H at the maximum.

.5 mm PbEq shield, 95.83% at 100 kV
Lead Shield measures 24" x 24"



SHS-LB-MB-2424

MOBILE LEAD PORTA-SHIELD

Includes durable, Nylon covered .5mmPbEq shield
Adjustable height from 36" to 56"

Available in three sizes:
24x24, 24x36, 24x48
Base: 24" Diameter

OVERHEAD ADJUSTABLE SHIELDS

CEILING MOUNTED OPTIONS



SHS-PTO-005

Ceiling Mounted Overhead Lead Acrylic Barrier With Lead Curtain and Light.

Lead Acrylic Window, 0.5mm PBEq
Window measures 24" x 16"
Lead Curtain, 0.50mm PbEq
Lead Curtain measures 24" x 14"



SHS-PTO-007

Ceiling Mounted Overhead Lead Acrylic Barrier With Lead Curtain
Lead Acrylic Window, 0.5mm PBEq
Window measures 24" x 16"
Lead Curtain, 0.50mm PbEq
Lead Curtain measures 24" x 14"



SHS-PTO-006

Ceiling Mounted Overhead Lead Acrylic Barrier With Torso Cutout
Lead Acrylic Window, 0.5mm PBEq
Window measures 32" x 24"
Torso cut-out to help fit most patients.



SHS-PTO-008

Ceiling Mounted Overhead Lead Acrylic Barrier With Torso Cutout and Light
Lead Acrylic Window, 0.5mm PBEq
Window measures 32" x 24"
Torso cut-out to help fit most patients. Examination light on second suspension arm.

ADDITIONAL PROTECTION

TABLE PADS, DRAPES, POSITIONERS
& MORE!




TABLE DRAPES





The single sided table drape uses velcro to adhere to the side of your imaging table. The double sided drape features a saddlebag design for a secure fit to any style imaging table. They provide protection against scatter radiation for your lower body and are specially designed to not obstruct positioning of the C-arm. Each drape attaches to the top by velcro allowing for easy cleaning and adjustability.

DRAP DIMENSIONS	MODEL # SINGLE SIDED	MODEL # DOUBLE SIDED
48" x 24"	SH-TS-4824	SH-TD-4824
48" x 28"	SH-TS-4828	SH-TD-4828
48" x 32"	SH-TS-4832	SH-TD-4832
48" x 36"	SH-TS-4836	SH-TD-4836
36" x 24"	SH-TS-3624	SH-TD-3624
36" x 28"	SH-TS-3628	SH-TD-3628
36" x 36"	SH-TS-3636	SH-TD-3636

SURGICAL TABLE SHIELDS

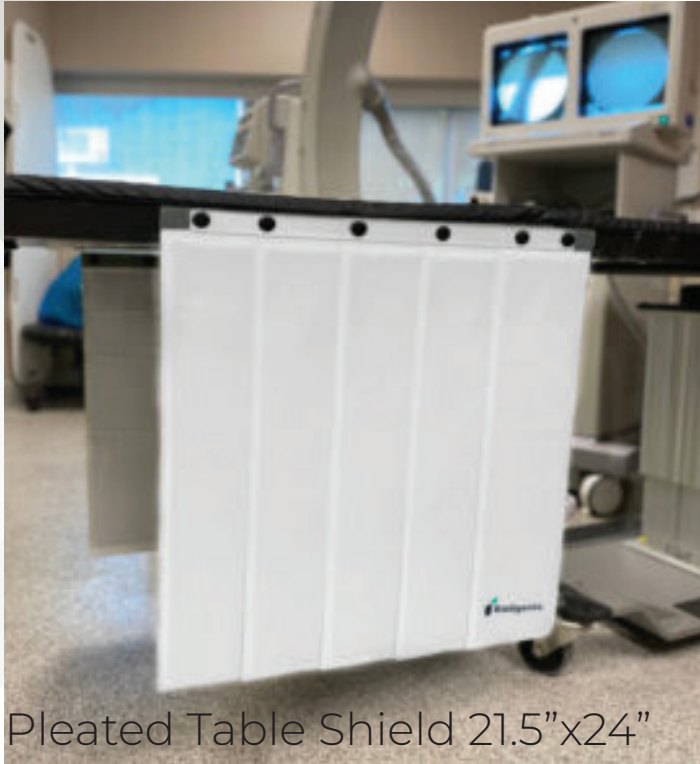


Pleated Table Shield 21.5"x24"

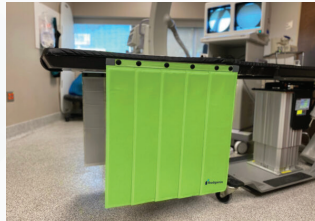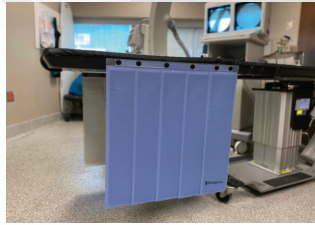

Radgenic shield adds an extra layer of protection with a lead core of 0.5mm Pb Eq. Available in white, blue, grey and green. Made from a polyurethane coated fabric which can be easily cleaned and disinfected with Quaternary or Isopropyl disinfectant wipes.

The Radiation Shielding Pleated Table Shield is non-porous, with a smooth surface, and fully sealed seams. It is also latex free, has no Velcro, and is soil and stain resistant.

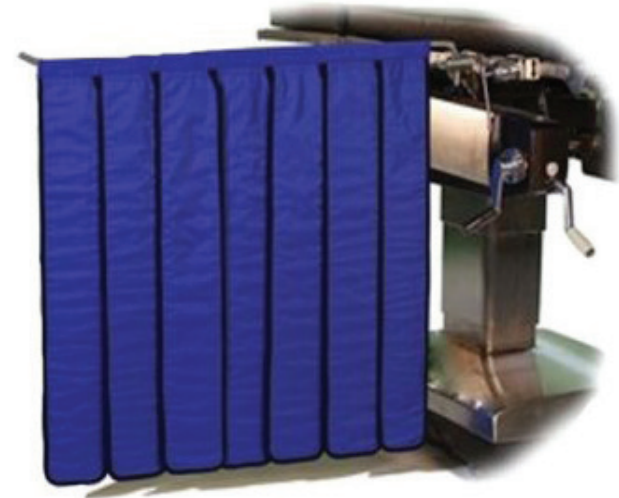

Lead Table Shield 36" x 36"


Lead table shield pivots a full 180 degrees to adjust to any angle you desire. Pleated "Pass Through" design to help facilitate any rotation of the C-ARM

The lead table shield attaches to the operating room table by inserting the stainless steel 5/8" rod into a standard side rail Clark socket. Clark socket not included.


Panel Size: 36" x 36"
.5mm lead equivalence

DISPOSABLE SURGICAL RADIATION PADS


Tested in accordance with ASTM and IEC standards, ScatterGuard drapes provide excellent protection against scatter radiation for fluorolab and operating suite personnel during interventional procedures. These drapes are placed upon the patient's body, so you can enjoy the benefits of added protection without any extra weight! They are sterile, disposable, conform to the patient's body and feature 3M adhesive backing to prevent slippage.

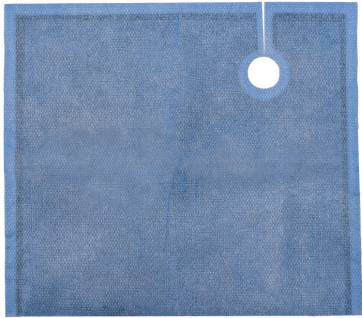


ANGIO DRAPE

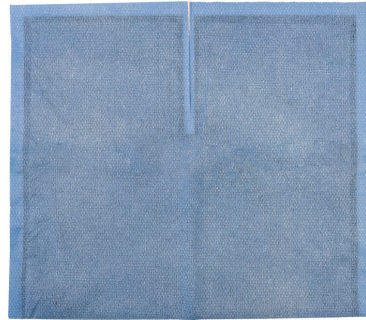


BILIARY SPLIT ACCESS

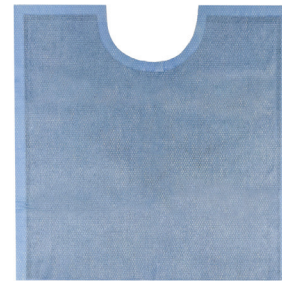


EP DRAPE LARGE

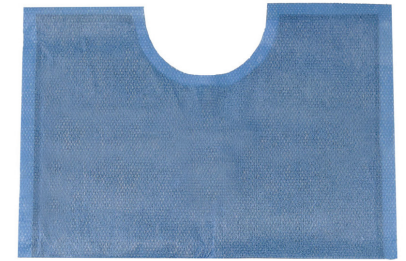


EP DRAPE SMALL

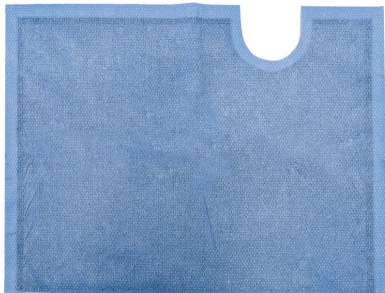

EP LEFT SUBCLAVIAN

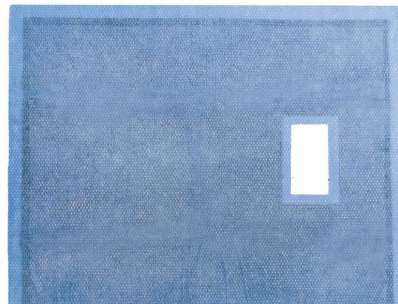


MULTI-PURPOSE WITH BOX

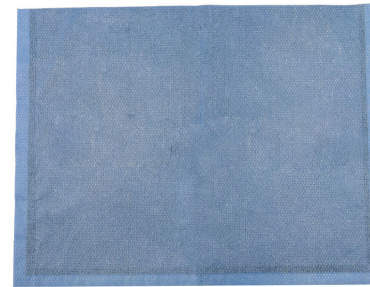


MULTI-PURPOSE



PERIPHERAL DRAPE

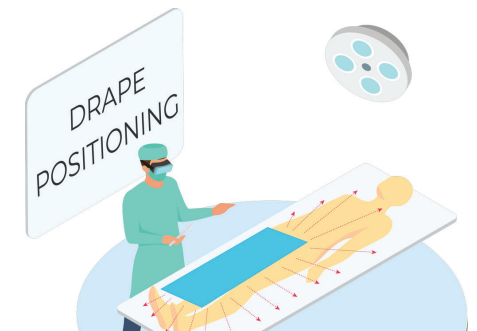

SURGICAL TABLE PADS

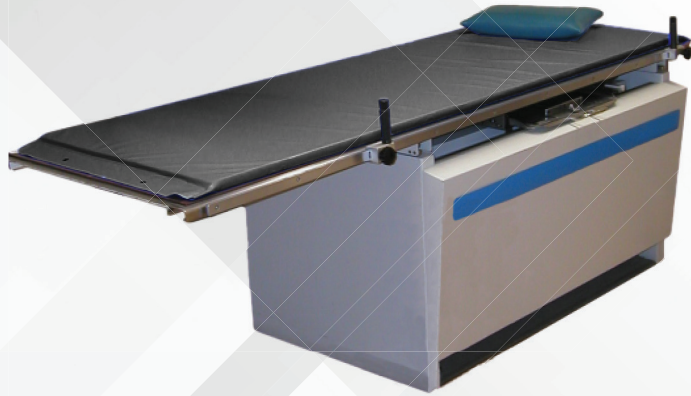

PREMIUM RADIOLUCENT X-RAY TABLE PADS



Table Pad Without Grommets



Table Pad with Grommets



ECONOMY RADIOLUCENT X-RAY TABLE PADS

- Comfortable and durable vinyl over a premium foam interior.
- Perfect for x-ray tables and stretchers.
- Increases patient comfort and stability.
- Available in standard sizes designed for virtually any application.
- Feature tapered ends with brass grommets for hanging storage.
- Radiolucent* and cleanable surface!

PREMIUM AND ECONOMY AVAILABLE SIZES:

Standard

72" x 23.25" x 1" Comfort Foam

72" x 23.25" x 2" Firm Foam

Standard Plus

80" x 30" x 1" Comfort Foam

80" x 30" x 2" Firm Foam

* Both Premium & Eco available with and without grommets

PREMIUM RADIOLUCENT X-RAY TABLE PADS



- Comfortable and durable fabric over a premium foam interior.
- Available with reinforced vinyl or 4-way stretch comfort fabric.
- Perfect for x-ray tables and stretchers.
- Increases patient comfort and stability.
- Premium boxed corners provide a superior finish and look.
- Available in standard sizes designed for virtually any application.
- Radiolucent* and cleanable surface!

AVAILABLE SIZES:

PREMIUM

72" x 23.25" x 1" Comfort Foam

72" x 23.25" x 2" Firm Foam

72" x 23.25" x 2" Memory Foam

PREMIUM PLUS

80" x 30" x 1" Comfort Foam

80" x 30" x 2" Firm Foam

80" x 30" x 2" Memory Foam

SURGICAL TABLE PAD REPLACEMENT




Tired of getting dings from Joint Commission for cracks, tears, rips or holes in your mattresses and surgical pads? Allow our team of trained experts to inspect your existing inventory and restore it or replace it!

Our high quality mattress pads and patient pads are made in the USA and offer premium quality with quick turnarounds. We have replacement pads for all styles of hospital equipment.

Stop paying TOP DOLLAR for name brand pads. We have Quality Mattress pads at affordable prices.

LEVELS OF PADDING AND COVERING

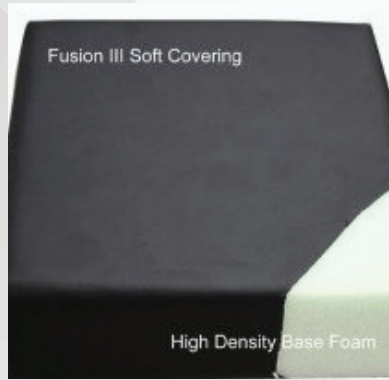

SHEERCARE SERIES LEVEL I

A great place to start! This surgical table pad is constructed of high density base foam surrounded by a Fusion III® ultra soft 4-way stretch covering. Our SheerCare Series offers an economical alternative without compromising the comfort and safety of your patients.

High Density Base Foam → Fusion III® Ultra Soft Covering



ULTRACARE SERIES LEVEL III

Our UltraCare Series surgical table pad sets take “Pressure Management” to the next level. These pad sets are constructed with high density base memory foam, then topped with a 1/2" layer of Visco memory foam and surrounded by a Fusion III® ultra soft 4-way stretch covering.

High Density Base Foam → Visco Memory Foam → Fusion III® Ultra Soft Covering



COMBO SERIES LEVEL III

Our Combo Series is our most comprehensively designed pressure management surgical table pad system. Our unique “tri-level” design is constructed of a high-density supportive base foam then topped with a 1/2" layer of doughy soft cradling Visco memory foam. A layer of Blue Diamond® Gel is placed just below the Fusion III® ultra soft 4-way stretch covering to ensure your patients the maximum protection available for pressure management.

High Density Base Foam → Visco Memory Foam → Blue Diamond Gel® → Fusion III® Ultra Soft Covering

COMMON SILHOUETTES

SKYTRON
STRYKER
AMSCO
PEDIGO
HILLROM
STERIS
SEIMENS
BERCHTOLD
MAQUET
AND MORE!

Our team of mattress experts
 are capable of making pads for
ALL types and brands of
 equipment.

AMSCO/STERIS
1080/2080



AMSCO/STERIS
3080



SKYTRON
6700,6701



SKYTRON
3500



MAQUET
ALPHASTAR



SHAMPAINE
4900/5100



BERCHTOLD
B710/B810



MAQUET 1130
SPLIT LEG



SHAMPAINE
1501 RADI-OP



MIDMARK
MONTREAL







STRETCHER PADS & ARMBOARDS



Measurement Guide

Indicate the style needed and provide the required measurements




Surgical Pads

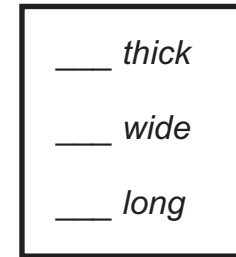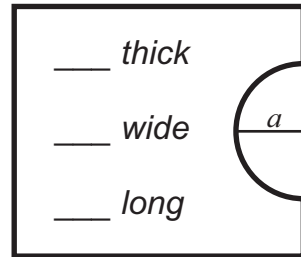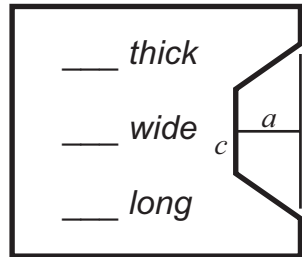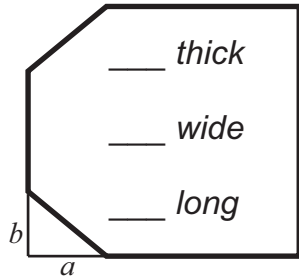


Select a Pad Style

- Economy
- Standard
- Premium

Stretcher Pads

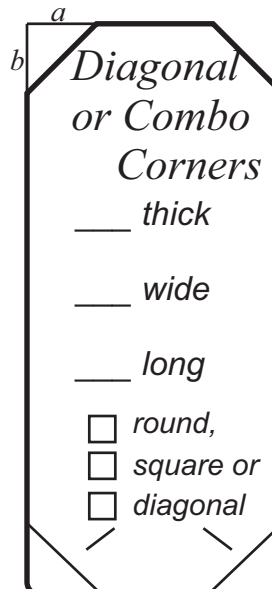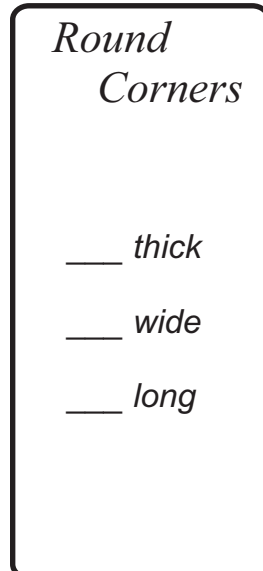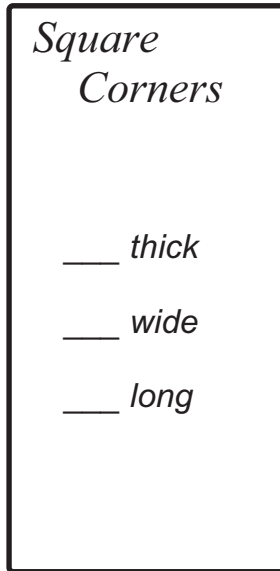

Armboard Pads

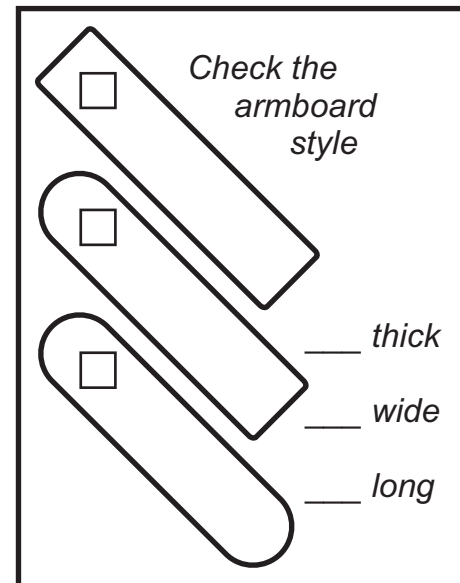


Corners & Cutouts

a _____
b _____
c _____

This is a 2" radius corner.

MARKERS

LEAD MARKERS FOR X-RAY

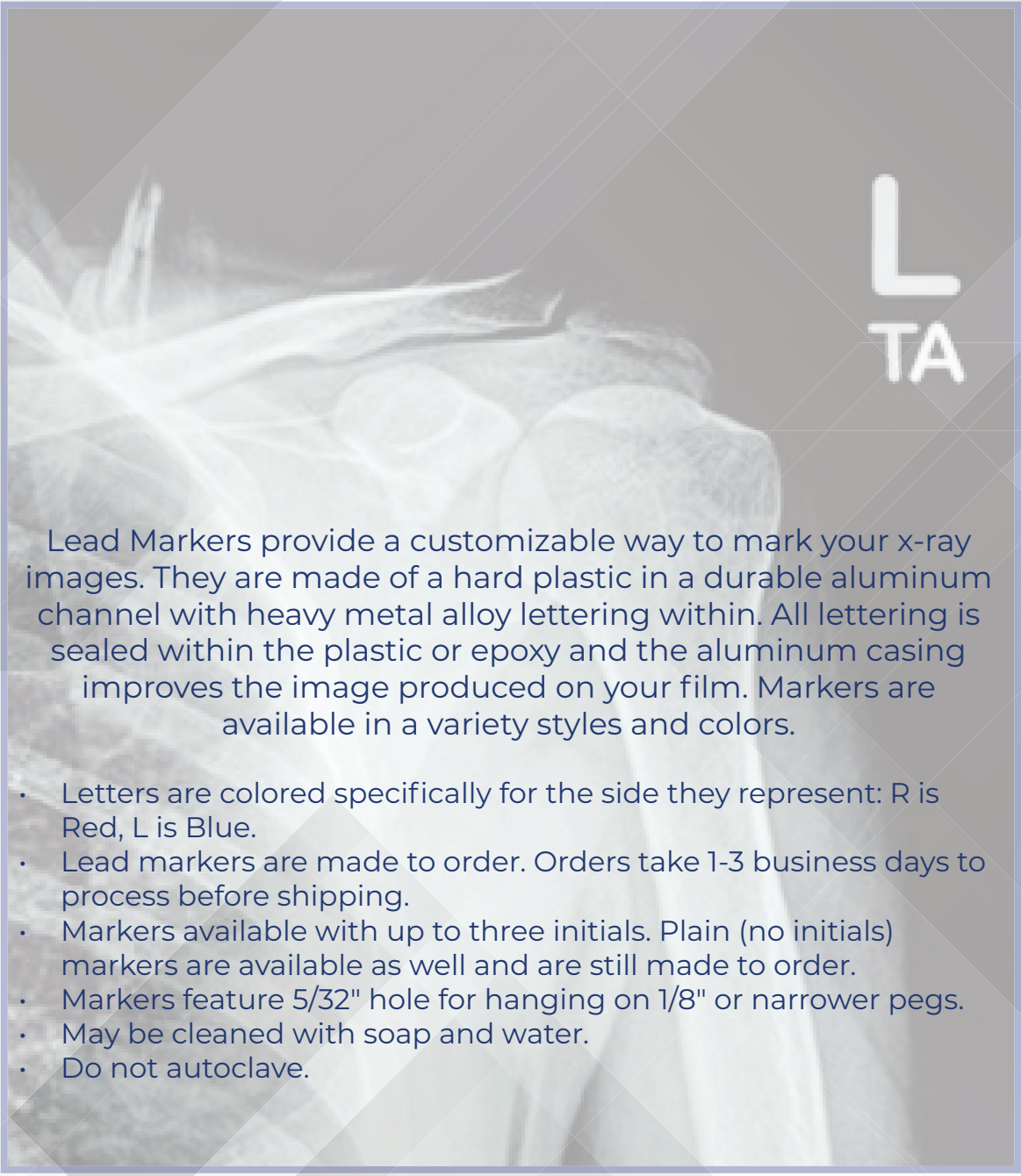An X-ray image of a knee joint. A white lead marker is placed on the right side of the image. The marker consists of a vertical bar with a horizontal bar extending to the right, forming an 'L' shape. Below the vertical bar, the letters 'TA' are printed. The X-ray shows the bones of the knee, including the femur, tibia, and patella.

Lead Markers provide a customizable way to mark your x-ray images. They are made of a hard plastic in a durable aluminum channel with heavy metal alloy lettering within. All lettering is sealed within the plastic or epoxy and the aluminum casing improves the image produced on your film. Markers are available in a variety styles and colors.

- Letters are colored specifically for the side they represent: R is Red, L is Blue.
- Lead markers are made to order. Orders take 1-3 business days to process before shipping.
- Markers available with up to three initials. Plain (no initials) markers are available as well and are still made to order.
- Markers feature 5/32" hole for hanging on 1/8" or narrower pegs.
- May be cleaned with soap and water.
- Do not autoclave.

LEAD MARKERS



SHS-TA-PPP11-V
Plastic (5/8") Vertical
Indicator Markers Set
(No Characters)



SHS-TA-PAP03-V-2
Aluminum (5/8") Vertical
Indicator Markers Set
(2 Lines - 3 Character Max)



SHS-TA-EAP003-2
Aluminum (5/8") Vertical
Marker Set
(2 Lines - 3 Character Max)



SHS-TA-EPP003
Plastic 5/8" Set
3 Character Max



SHS-TA-PPP13-V
Plastic (5/8") Vertical
Indicator Markers Set
(3 Character Max)



SHS-TA-EPP001
Compact Plastic 5/8" Set
(No Characters)



SHS-TA-EAOP01
Aluminum Oncology
(3/4") Vertical Marker Set
(No Characters)



SHS-TA-EPP013-2
Plastic 5/8" Set
(2 Line - 3 Character Max)



SHS-TA-EPP013
Plastic (1") Vertical Marker
Set (3 Character Max)



SHS-TAEAP003
Aluminum 1/2" Set
3 Character Max




SHS-TA-EAP013
Aluminum (1") Vertical
Marker Set (With 3
Character Max)



SHS-TA-EAP003-625
Aluminum 5/8" Set
3 Character Max



SHS-TA-EAP011
Aluminum (1") Vertical
Marker Set (No
Characters)

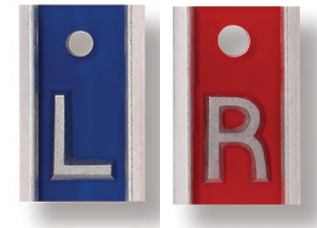

SHS-TA-EAP001
Aluminum (5/8")
Vertical Marker Set (No
Characters)

LEAD MARKERS CONTINUED



SHS-TA-PAP03-V
Aluminum (5/8") Vertical
Position Indicator Marker
Set (3 Character Max)



SHS-TA-PPP01-V
Aluminum (5/8") Vertical
Position Indicator Marker
Set (No Characters)



SHS-TA-PPP03-V
Plastic L/R Compact
Vertical 7/8" Indicator Set
(3 character max)



SHS-TA-PPP01-V
Plastic L/R Compact 7/8"
Set W/ Indicator



SHS-LM-001-A



SHS-LM-007




SHS-LM-004




SHS-LM-003



SHS-LM-006

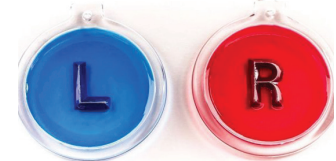

SHS-LM-009

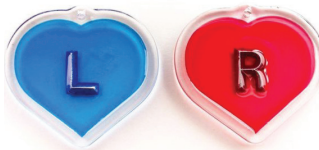

SHS-LM-010



SHS-LM-012



SHS-LM-005



SHS-LM-008




SHS-LM-011

LEAD MARKERS CONTINUED



SHS-TA-EPA1000
Plastic Arrow Multi-
Directional Marker (No
Characters)

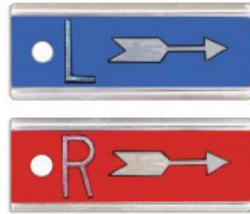

SHS-TA-EAP041
Aluminum (1/2")
Horizontal Arrow
Marker Set (No
Characters)



PLASTIC WORD
MARKERS



ALUMINUM WORD
MARKERS




SHS-TA-EAP021
Plastic (1/2") Horizontal
Marker Set (No
Characters)



SHS-TA-EA1000
Aluminum Arrow Multi-
Directional Marker (No
Characters)



SHS-TA-EPP032
Plastic (1/2") Horizontal
Arrow Marker Set (With 3
Character Max)

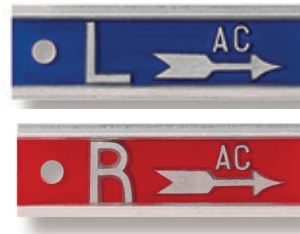

SHS-TA-EAP043
Aluminum (1/2")
Horizontal Arrow Marker
Set (With 3 Character
Max)



SHS-TA-EAP032
Aluminum (1/2")
Horizontal Marker Set
(With 3 Character Max)



SHS-TA-EAP021
Aluminum (1/2")
Horizontal Marker Set (No
Characters)



SHS-TA-PAP03-H
Aluminum (1/2")
Horizontal Position
Indicator Markers
(With 2 Character Max)



Contact your local sales representative for customized lead markers. These are available in a variety of configurations, styles, colors and materials. For more information call (800) 674-3871

MARKER COLOR OPTIONS (SOLID)



RED
(RIGHT DEFAULT)



BRILLIANT BLUE
(LEFT DEFAULT)



DARK BLUE



COBALT BLUE



ICE BLUE



TURQUOISE



DARK GREEN



YELLOW
GREEN




LIME TREE
GREEN



BRIMSTONE
YELLOW



ORANGE



BROWN



VIOLET



PURPLE RED



BURGUNDY



DARK RED



PINK



SOFT PINK



CREAM



GOLD
METALLIC




SILVER
METALLIC



GREY



BLACK




WHITE


MARKER COLOR OPTIONS (GLITTER)



RED
(RIGHT DEFAULT)

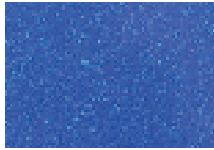

BLUE
(LEFT DEFAULT)

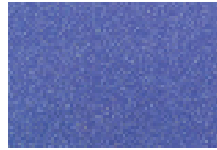


LIGHT BLUE



FLUORESCENT
BLUE



TIFFANY BLUE

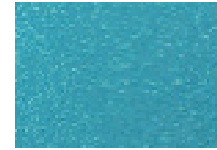

MINT



SEA FOAM



GREEN




FLUORESCENT
GREEN



ORANGE




FLUORESCENT
ORANGE

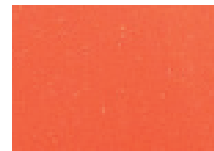

CORAL




PURPLE

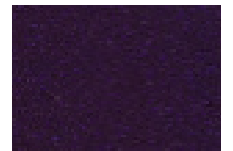

DARK
AMETHYST




DARK RED

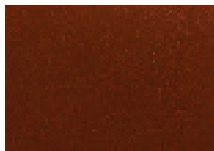


CINNAMON



FLUORESCENT
PINK

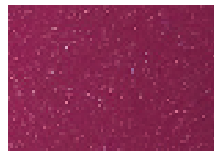


MAGENTA

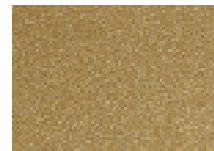


GOLD




ROSY GOLD

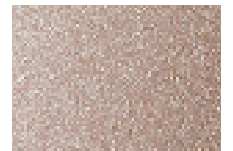


CHAMPAGNE

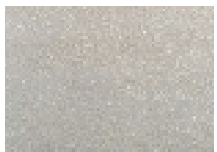

SILVER

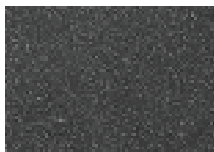

DARK GRAY

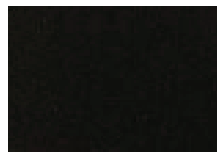

BLACK



WHITE


SPONGES & POSITIONING




NO MATTER WHAT PROCEDURE, OUR EXTENSIVE LINE OF PATIENT POSITIONING AIDS AND SPONGES HAS YOU COVERED. CREATE A CUSTOM BUNDLE, ORDER REPLACEMENTS OR CHOOSE FROM ONE OF OUR BUNDLES ALREADY FORMULATED FOR THE OPERATING ROOM NEEDS.

- CLOSED CELL SPONGES
- COATED SPONGES
- NON-COATED SPONGES
- VINYL COVERED SPONGES
- MRI POSITIONING AIDS
- SPECIALTY POSITIONING AIDS
- PLATFORMS & POSITIONING STOOLS
- SANDBAGS & SANDSACKS

CLOSED CELL SPONGES

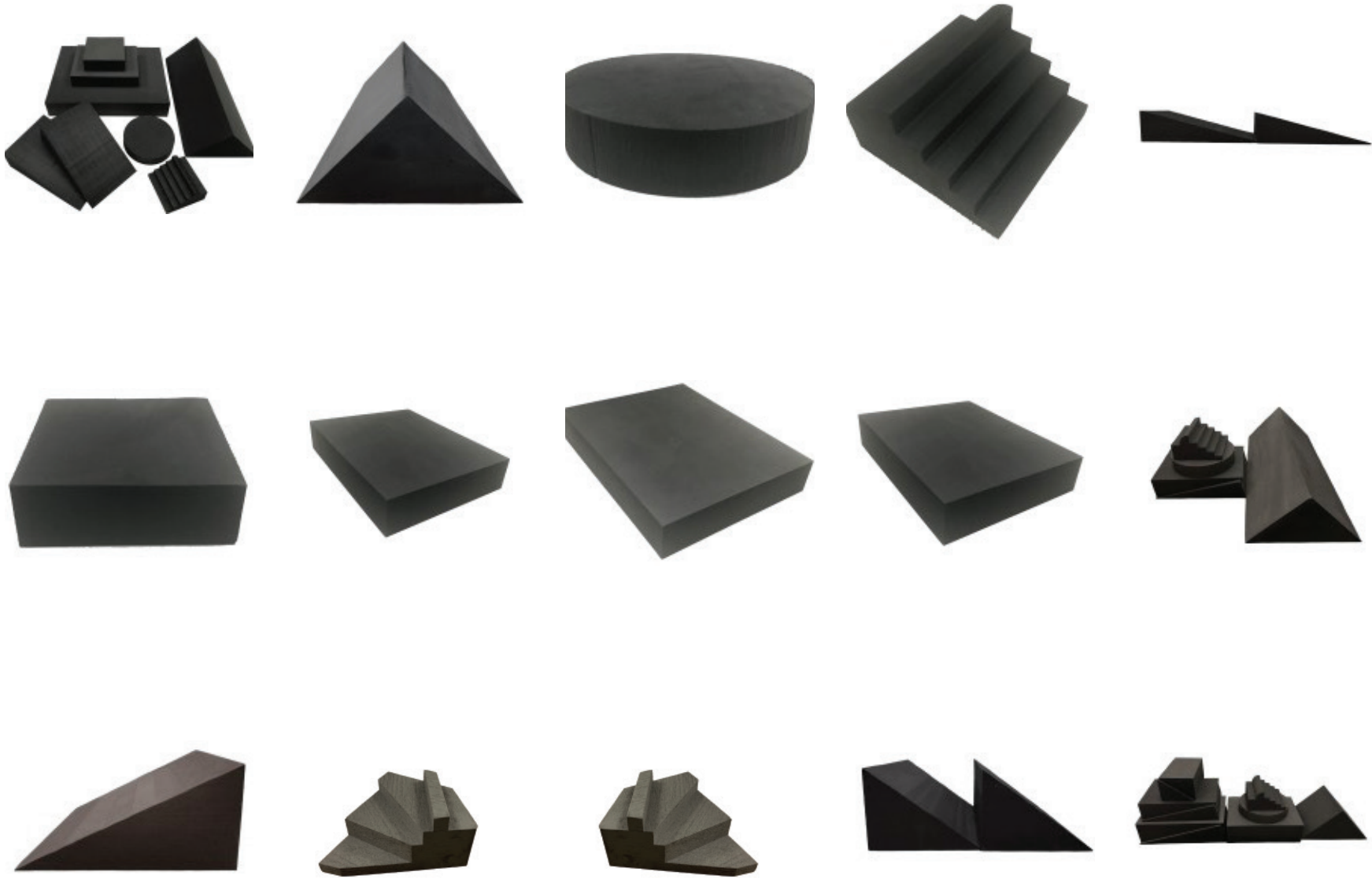

CLOSED CELL SPONGES - CONTINUED

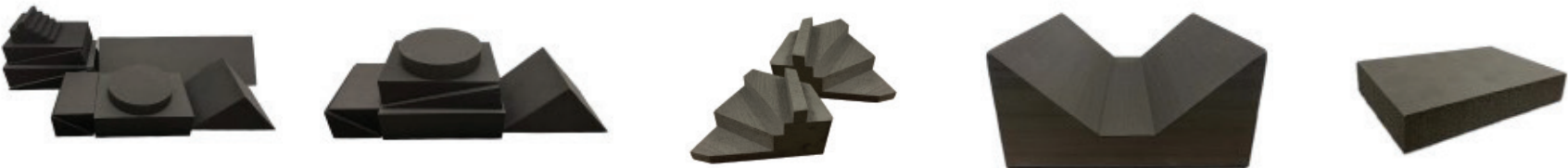

CLOSED CELL - CONTINUED



Closed Cell Sponges are made from a lightweight, radiolucent, non-porous foam. The perfect solution for any positioning need! These positioning sponges are entirely fluid resistant, easy to clean and extremely durable! Featured in a non-coated, premium charcoal grey color.

- 100% Radiolucent
- Artifact free imaging
- Cleanable, fluid resistant surface
- Lightweight and comfortable
- The number one option in positioning sponges!

VINYL COVERED SPONGES



VINYL COVERED SPONGES - CONTINUED




Vinyl Covered Sponges are made from a high density, lightweight foam and covered in our royal blue, heavy gauge, reinforced vinyl. These premium positioning sponges help create extreme patient comfort during long procedures! Featuring a non-porous, cleanable surface, our durable vinyl sponges also include completely sewn seams to ensure complete protection and sanitation!


- Helps ensure extended patient comfort
- Extremely durable and lightweight
- Cleanable and fluid resistant premium vinyl surface
- Completely sewn seams ensure complete sanitation!

Vinyl covered sponges are radiolucent but seams will show under x-ray.

COATED SPONGES - NON-STEALTH

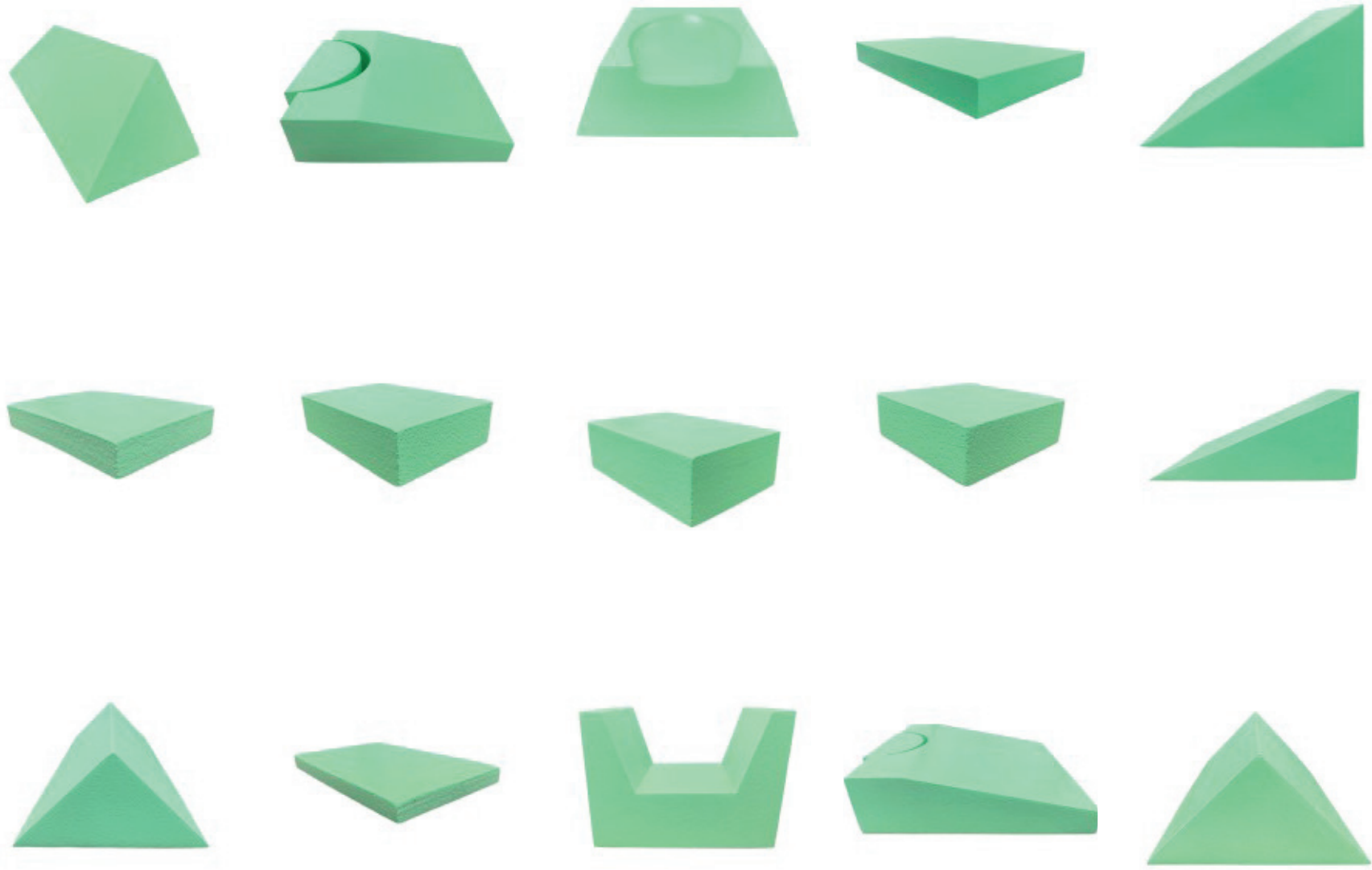

COATED SPONGES - NON-STEALTH

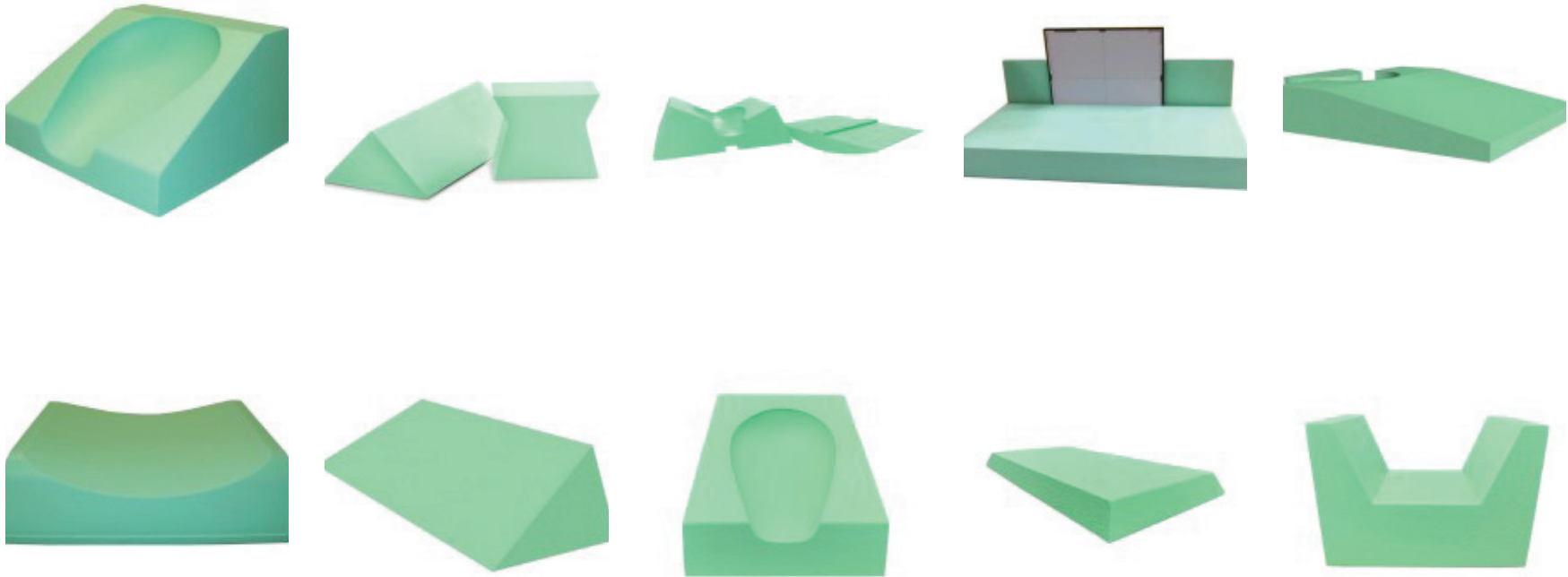


PEDIATRIC

Coated Sponges feature a proprietary medical grade, fluid-resistant, cleanable coating over a standard (core) foam interior. The coating applied protects the sponge from the absorption of body fluids. Minimal shadowing with low attenuation. Light green color. Coated sponges are radiolucent.

Non-Stealth Sponges are formed with traditional angles and do not feature our proprietary “flared” angles to offer significantly reduced artifacts during viewing. Non-stealth (or non-flared) sides may show artifacts when imaging, specifically when coating is applied. A variety of bundles are available. From pediatric to bariatric patients, we offer the largest selection of positioners on the market!


Anytime a coating is applied to sponges it will show under x-ray.

COATED SPONGES - STEALTH

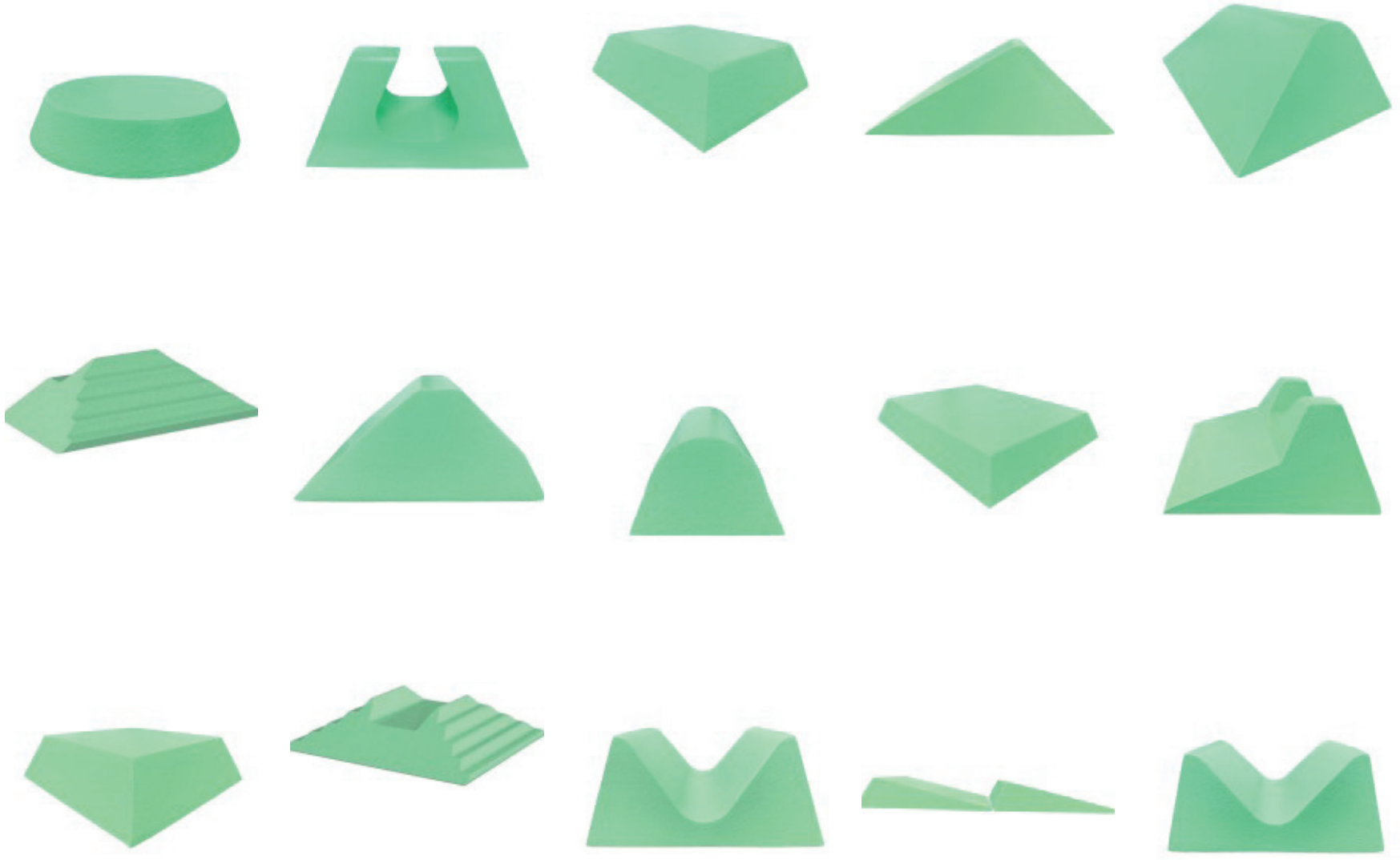


PEDIATRIC

COATED SPONGES - STEALTH

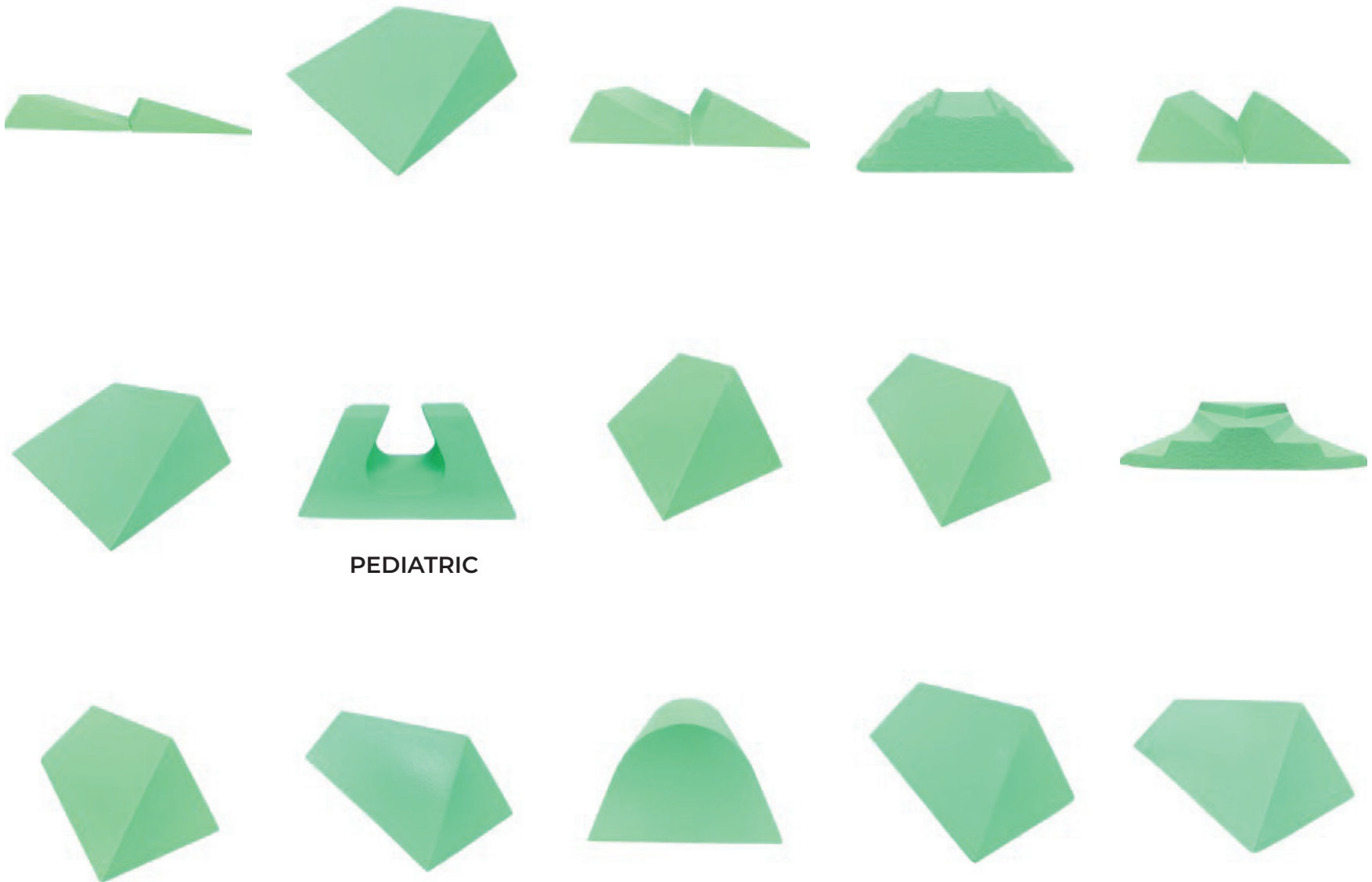

PEDIATRIC

COATED SPONGES - STEALTH

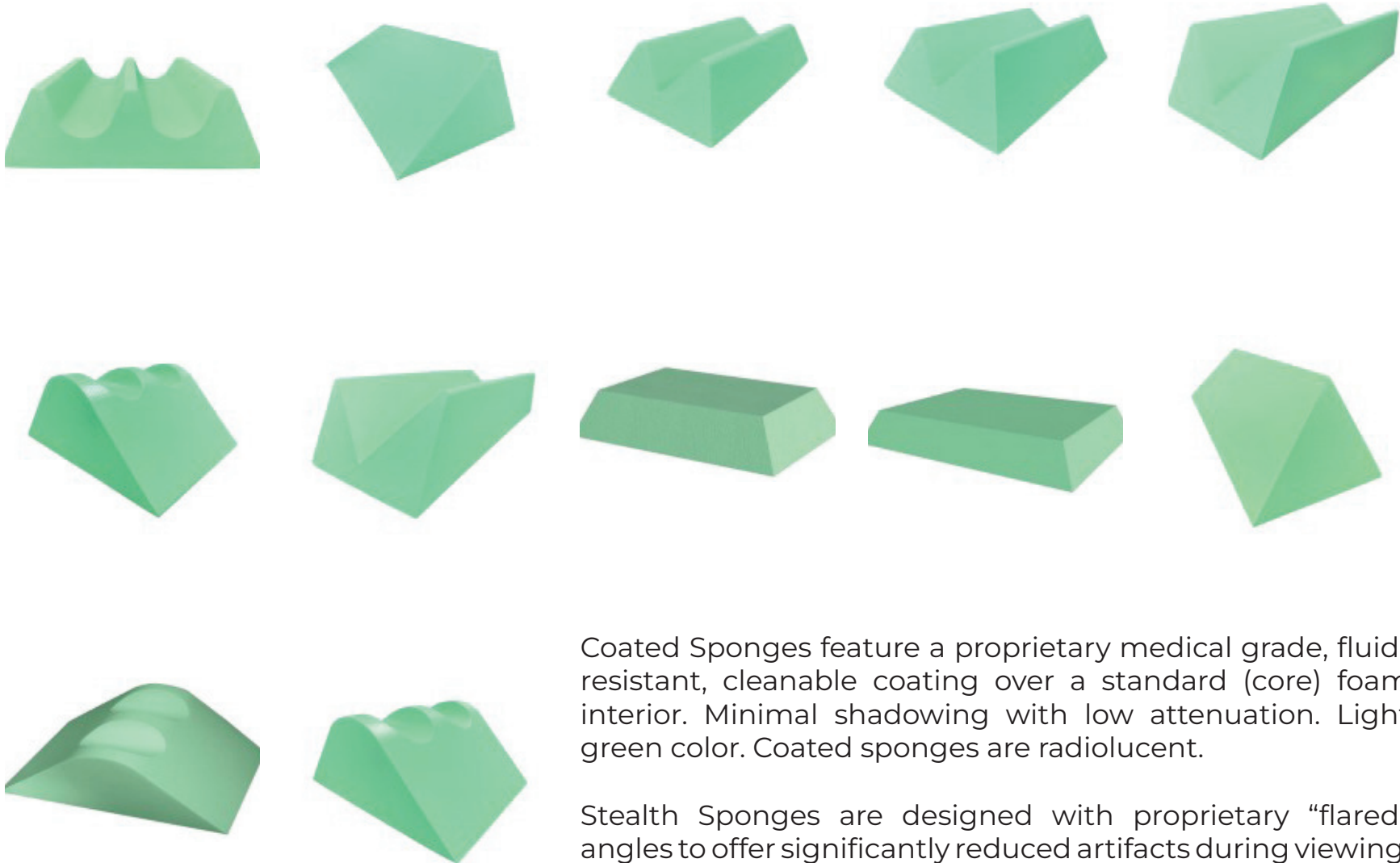

Coated Sponges feature a proprietary medical grade, fluid-resistant, cleanable coating over a standard (core) foam interior. Minimal shadowing with low attenuation. Light green color. Coated sponges are radiolucent.


Stealth Sponges are designed with proprietary “flared” angles to offer significantly reduced artifacts during viewing. Flared sides will increase overall size.

MRI NYLON COVERED SPONGES



MRI Positioners meet all requirements for MRI compatibility and safety. Approved for use up to 3 Tesla. The unique viscoelastic material is specifically designed for patient comfort and durability. The coated nylon is flexible, fluid & stain resistant and features permanent anti-microbial properties!

- Approved for use in MRI (up to 3 Tesla)
- Extreme comfort for patients
- Lightweight and durable
- Permanent anti-microbial properties
- Breathable waterproof polymer fabric
- Stain resistant and latex free

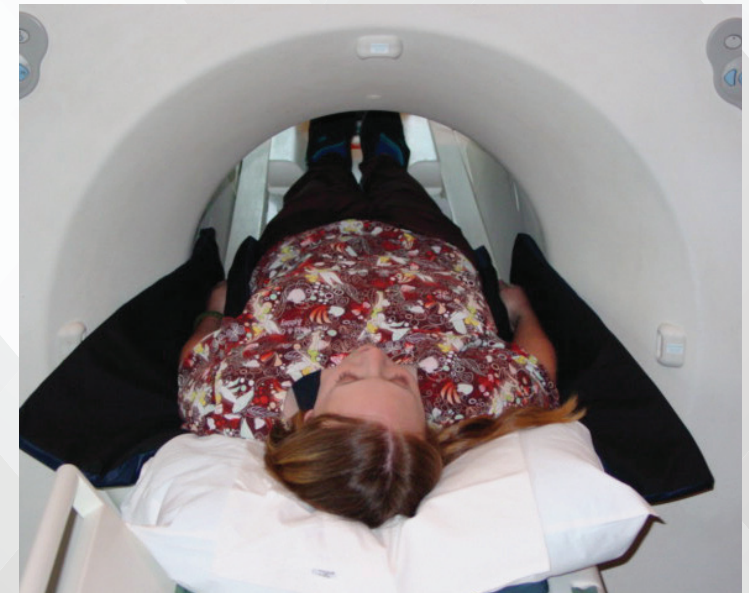

CUSTOM BUNDLES AVAILABLE



BUNDLES ARE AVAILABLE IN CLOSED CELL, COATED AND NON-COATED - STEALTH AND NON-STEALTH, MRI VINYL COVERED. ALL BUNDLES NOT SHOWN. CALL US TODAY FOR MORE OPTIONS AND INFORMATION

800-674-7831

BLUE DIAMOND GEL POSITIONERS



- Radiolucent
- Antimicrobial & Antibacterial
- Reusable, Repairable, and Easily Cleaned
- Latex, Silicone, Plasticizer Free
- Contains No Water, Will not Leak, Ooze, or Evaporate
- Redistributes Weight & Provides Pressure Reducing Attributes
- 2 Year Warranty

SANDBAG POSITIONERS + BUNDLES



Our premium, heavy-gauge vinyl sandbags are the ideal option for quick and comfortable patient positioning. These sandbags are the perfect solution for nearly every positioning need such as: veterinary use, therapeutic and rehabilitation applications, surgical after-care, and general patient positioning. Each sandbag features a heat sealed inner sleeve, covered by a completely sewn and reinforced vinyl for long lasting durability.

Available with either a standard comfort loop handle or an upgraded hook & loop, AC joint wrist strap. Choice of Black or Navy Vinyl. Available in 5, 7 and 10 pound weight.

Approved for use in MRI (up to 3 Tesla)

- Extremely durable, reinforced vinyl.
- Heat sealed inner sleeves prevent leakage.
- Completely sewn, tear resistant seams.
- Designed to meet nearly any positioning need.
- Available with either a standard comfort loop handle or an upgraded AC joint wrist strap.

FLEXIBLE VINYL SANDSOCKS + BUNDLES



SINGLE WEIGHTS

SS-01
SS-03
SS-04
SS-05
SS-07
SS-08
SS-10
SS-11
SS-12
SS-15



GENERAL BUNDLE 1

SS-01 - 1
SS-03 - 1
SS-05 - 1
SS-10 - 1
SS-12 - 1



GENERAL BUNDLE 2

SS-01 - 2
SS-03 - 2
SS-05 - 4
SS-11 - 2
SS-12 - 4
SS-15 - 2



GENERAL BUNDLE 3


SS-01 - 2
SS-03 - 2
SS-05 - 4
SS-11 - 2
SS-12 - 4
SS-15 - 2

We offer superior, flexible vinyl sandsocks are the ideal option for quick and comfortable patient positioning. These flexible sandsocks are the perfect solution for hard to hold or irregular anatomy. Designed for nearly every positioning need such as: veterinary use, therapeutic and rehabilitation applications, surgical after-care, and general patient positioning. Each tan colored vinyl sandsock features reinforced double seams for extended durability. Available in Tan (Camel) or Black vinyl and in weight increments of: 1, 3, 4, 5, 7, 8, 10, 11, 12, 15 pounds.

KEY FEATURES


Reinforced double seams for extended durability
Extremely flexible, tan colored vinyl covering
Lightweight and easy to use
Designed to meet nearly any positioning need

SPECIALTY PATIENT RESTRAINT

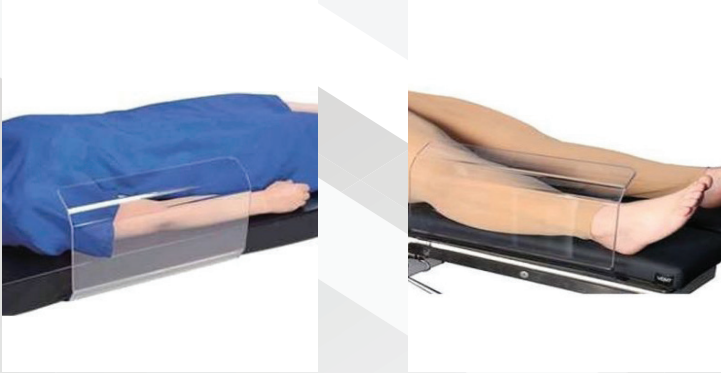

Toboggan Arm Boards are easy-to-clean radiolucent shields made of a crystal clear, extremely strong, and radiolucent polycarbonate material. They protect the patient's arms/legs during interventional procedures. Slide the toboggan under the pad and the patient's weight holds it in place. Contoured top offers stabilizing support for clinician's arms while shielding against sharp edges during procedures. MRI safe. Available in: SM: 10"L x 6"H x 10-1/4"D
MED: 15-3/4"L x 8"H x 10-1/4"D and LG: 15-3/4"L x 10-1/2"H x 10-1/4"

The Radiolucent Papoose Board is a MRI-Safe tool for infants to adults. Unlike alternatives that use metal, which creates artifacts in the scans, this frame is uniquely metal-free. As a result, you will receive clear scans during the diagnostic imaging process. This complete set includes the board, a padded three-tier flap set, a padded head strap, as well as 4 arm/wrist straps. Available in SMALL, REGULAR, LARGE and XL. *regular size shown

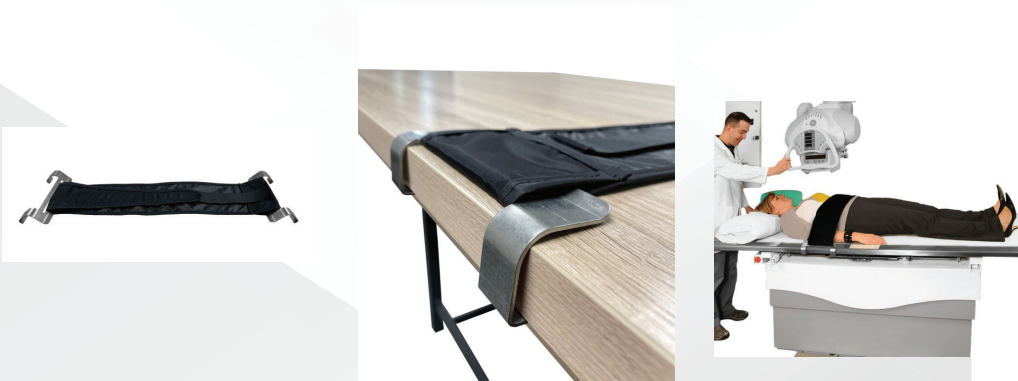

The Universal Patient Restraint Strap is the perfect solution for securing a patient during imaging. This 52" x 6" universal nylon strap is specifically designed to fit nearly any imaging table. Each strap fits quickly, firmly and comfortably over patients during imaging. Made of a high quality, durable black nylon with a convenient and easy to adjust hook & loop closure. The included 1.33" stainless steel bracket clamps for increased durability and strength.

PATIENT BOARDS



Specs:	5021	5051
Length	23"	30"
Width	12.5"	12.5"
Thick	1/4"	1/4"
Product Weight	2lbs	3lbs
Weight Cap.	300lbs	300lbs

Plastic Transfer Board - white - 2 sizes available



- Ultra-smooth sheet
- Anti-static and radiolucent
- High-density polyethylene
- Bridges gap between transfer points
- Six built-in handles
- Dimensions: 31" L x 23" W (tapers to 16" W) x 1/8" thick
- Weight Capacity: 250lbs

Slide Board Patient Transfer Sheet - white



- Semi-rigid polyethylene board is covered with an anti-static coating
- Ten handles for ease of positioning.
- Radiolucent
- Rounded edges
- MRI Safe
- Weight capacity: 250 lbs
- 18" W x 72" H x 3/16" thick
- Item Weight: 9 lbs

Static Treated Aqua Transfer Board - 18"



- Semi-rigid polyethylene board is covered with an anti-static coating
- Twelve handles for ease of positioning.
- Radiolucent
- Rounded edges
- MRI Safe
- Weight capacity: 250 lbs
- 22" W x 72" H x 3/16" thick
- Item Weight: 9 lbs


Static Treated Aqua Transfer Board - 22"

PATIENT BOARDS




- Semi-rigid polyethylene board is covered with an anti-static coating
- Twelve handles for ease of positioning.
- Radiolucent
- Rounded edges
- MRI Safe
- Weight capacity: 250 lbs
- 22" W x 72" H x 3/16" thick
- Item Weight: 9 lbs

Static Treated White Transfer Board

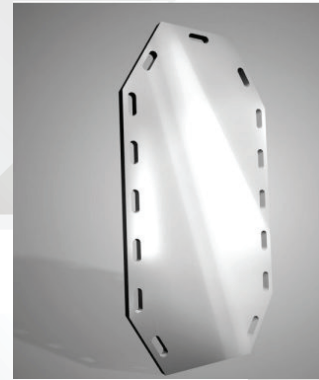


- Semi-flexible moving platform.
- Handle surrounds for ease of positioning.
- Radiolucent
- 1/4" thick moisture-proof high density polyethylene
- MRI Safe
- Weight capacity: 250 lbs
- 72" L x 32" W x 1/4" Thick

Bariatric Patient Transfer Board - white

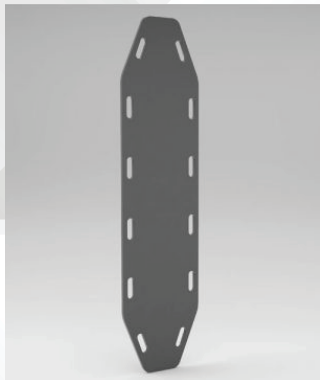

- Semi-rigid polyethylene board is covered with an anti-static coating
- Ten handles for ease of positioning.
- Radiolucent
- Rounded edges
- MRI Safe
- Weight capacity: 250 lbs
- 18" W x 72" H x 3/16" thick
- Item Weight: 9 lbs

Anti-Static Grey Transfer Board - 18"

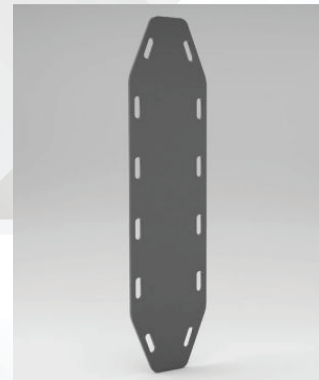

- Semi-rigid polyethylene board is covered with an anti-static coating
- Twelve handles for ease of positioning.
- Radiolucent
- Rounded edges
- MRI Safe
- Weight capacity: 250 lbs
- 22" W x 72" H x 3/16" thick
- Item Weight: 9 lbs

Anti-Static Grey Transfer Board - 22"

PATIENT BOARDS



- Transfers up to 600 lbs
- Reduces caretaker back strain
- Extra-wide design better supports excess tissue. It even supports patient's head.
- Fourteen handles
- Flexible plastic.
- 38"W x 72"L x 3/16" thick
- Weight: 16lbs

Bariatric Patient Shifter - Aqua



- Transfers up to 600 lbs
- Anti-Static Coating
- Extra-wide design better supports excess tissue. It even supports patient's head.
- Fourteen handles
- Flexible plastic.
- 38"W x 72"L x 3/16" thick
- Weight: 16lbs

Anti-Static Bariatric Patient Transfer Board



Product Specifications:

Specification	Dimensions
Width (center)	10.5"
Width (end)	12"
Length	28"
Product Weight	2.6 lbs
Weight Capacity	250 lbs.

Curved Banana Transfer Board



- Wall Mountable Rack adjusts to hold transfer board vertically or horizontally.
- One board can be stored per rack
- 18" L x 2.5" H x 1" D

Patient Transfer Board Storage Rack



STREAMLITE LEAD

Protection You Deserve

800-674-7831

STREAMLITELEAD.COM